



Recombinant Cluster of Differentiation 79B (CD79B)

Catalog No.: TP01534

50µg

Sequence Information

Species: Human

Gene ID:974

Swiss Prot:P40259

Synonyms:Igb; B29; Immunoglobulin-associated
B29 protein; B-cell antigen receptor
complex-associated protein beta chain

Residues:Asn37~Pro226

NPKGSACSRIWQSPRFIARKRGFTVKMH CYMNSASGNVSWLWKQEMDENPQQLK
LEKGRMEESQNESLATLTIQIRGFEDNGIYFCQKCNNTSEVYQCGGTEL RVMG
FSTLAQLKQRNTLKDGIIMIQTLLIILFIIVPIFLLDKDDSKAGMEEDHTYEG
LDIDQTATYEDIVTLRTGEVKWSVGEHP

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His Tag

Subcellular Location: Cell membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.4

Predicted Molecular Mass: 25.6kDa

Accurate Molecular Mass: 26kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

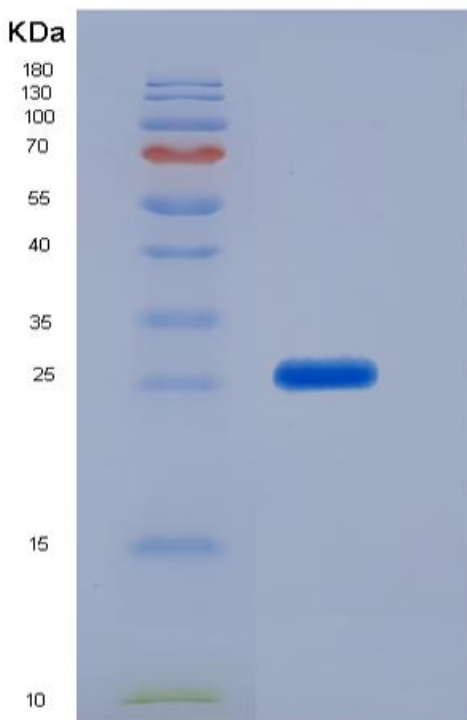


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.