



Recombinant Connective Tissue Growth Factor (CTGF)

Catalog No.: TP01765

50µg

Sequence Information

Species: Human

Gene ID:1490

Swiss Prot:P29279

Synonyms:CCN2; HCS24; IGFBP8; NOV2;

Insulin-Like Growth Factor-Binding

Protein 8; Hypertrophic

Chondrocyte-Specific Protein 24; CCN

Family Member 2

Residues:Gln27~Ala349

QNCSGPRCPDEPAPRCPAGVSLVLDGCGCCRVCAQLGELCTERDPCDPHKGL

FCHFSPANRKIGVCTAKDGAPCIFGGTVYRSGESFQSSCKYQCTCLDGAVGCM

PLCSMDVRLPSPDCPFPRRVKLPGKCCEEWCDEPKDQTVVGPALAAAYRLEDTF

GPDPMTIRANCLVQTTEWSACSKTCGMGISTRVTNDNASCRLEKQSRLCMVRPC

EADLEENIKKGKKCIRTPKISKPIKFELSGCTSMKTYRAKFCGVCTDGRCTPH

RTTTLPVEFKCPDGEVMKKNMFIKTCACHYNCPGDNDIFESLYYRKMVGDMA

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His and GST Tag

Subcellular Location: Secreted, Extracellular matrix.

Purity: >80%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.9

Predicted Molecular Mass: 65.5kDa

Accurate Molecular Mass: 65kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

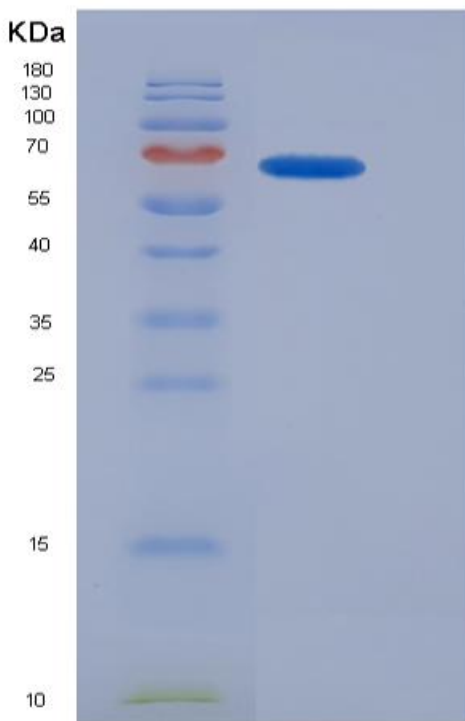


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.