

## Recombinant Tissue Factor (TF)

Catalog No.: TP08939

50µg

### Sequence Information

**Species:** Mouse

**Gene ID:**14066

**Swiss Prot:**P20352

**Synonyms:**CD142; FIII; F3; TFA;

Thromboplastin; Coagulation Factor III

**Residues:**Gly30~Glu251

GIPEKAFNLTWISTDFKTIIEWQPKPTNYTYTVQISDRSRNWKNKCFSTTDTEC  
DLTDEIVKDVWAYEAKVLSVPRRNSVHGDGDLVIHGEEPPFTNAPKFLPYRD  
TNLGGPVIQQFEQDGRKLNVVVKDSLTLVRKNGTFLTLRQVFGKDLGYIITYRK  
GSSTGKKTNITNTNEFSIDVEEGVSYCFFVQAMIFSRKTNQNSPGSSTVCTEQW  
KSFLGE

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:**N-terminal His Tag

**Subcellular Location:** Membrane.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:**PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

**Original Concentration:** 150µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 7.7

**Predicted Molecular Mass:** 26.6kDa

**Accurate Molecular Mass:** 31kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ IDENTIFICATION ]**

TGGATTTCGAGGAMGGSTTTAATTTACTTGGATACACTGATTTTCMAGICATTTTGGAGTGGGAGCCAAACCCACACACTATTCCTTCCTGTACGATTAAGTATGATCTAGAACTGGAAAGAGTCTCTCTGACGAGACCCAGTGGGACCTACAGAGGATGCTGAGGACTATGAGC  
 GIPEKAFNLTIVISTDFKTIILEWQPKPTNVTYTVQISDRSRNWKKCFSTTDTTECDLITDEIVKDVTVAYEA

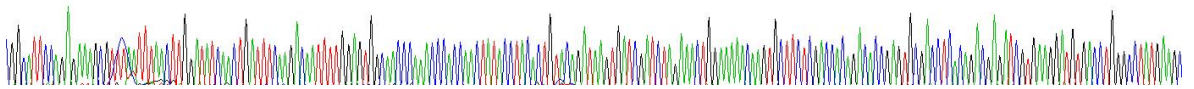


Figure 1. Gene Sequencing (Extract)

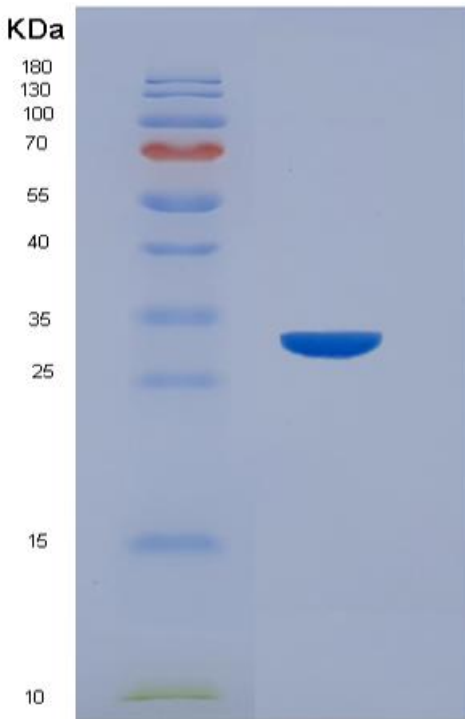


Figure 2. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any



issue if the kit was used in clinical diagnostic or any other procedures.