



# Recombinant General Transcription Factor IIA, Polypeptide 1 (GTF2A1)

Catalog No.: TP02367 50µg

## Sequence Information

**Species:** Human

**Gene ID:**2957

**Swiss Prot:**P52655

**Synonyms:**GTF2-A1; TF2A1; TFIIA;

Transcription Initiation Factor IIA

Subunit 1; Transcription initiation factor

TFIIA 42 kDa subunit

**Residues:**Ala2~Asp274

ANSANTNTVPKLYRSVIEDVINDVRDIFLDDGVDEQVLMELKTLWENKLMQSRA  
VDGFHSEEQQLLLQVQQQHQPQQQHHHHHHHQAQPQQTVPQQAQTQQVLIPA  
SQQATAPQVIVPDSKLIQHNMNASNMSAAATAATLALPAGVTPVQQILTNSGQLL  
QVVRANGAQYIFQPQQSVVLQQQVIPQMMPGGVQAPVIQQVLAPLPGGISPQT  
GVIIQPQQILFTGNKTQVIPTTVAAPTPAQAITATGQQQPQAQPAQTQAPLVL  
QVD

## Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:**Two N-terminal Tags, His-tag and GST-tag

**Subcellular Location:** Nucleus.

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:**PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:** 59.6kDa

**Accurate Molecular Mass:** 59kDa as determined by SDS-PAGE reducing conditions.

## **[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## **[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

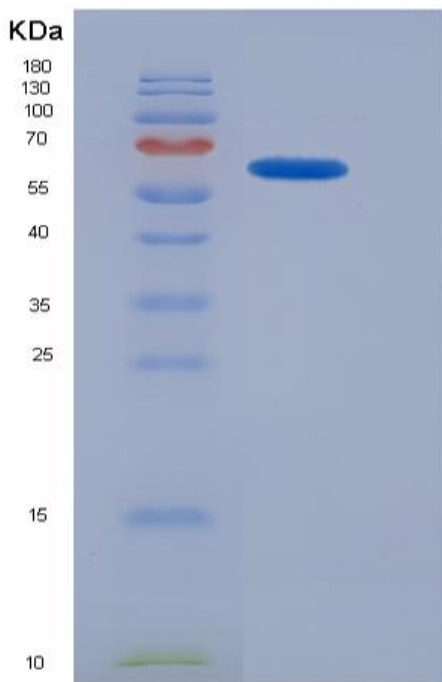


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.