

## Recombinant Gamma-Glutamyltransferase 1 (gGT1)

Catalog No.: TP08679

50µg

### Sequence Information

**Species:** Mouse

**Gene ID:**14598

**Swiss Prot:**Q60928

**Synonyms:**CD224; G-GT1; GTG;

Gamma-glutamyltranspeptidase 1;

Glutathione hydrolase 1; Leukotriene-C4

hydrolase;

**Residues:**Asn296~Pro537

NFSPKSVATPEQKALTYHRIVEAFRFAYAKRTMLGDPKFVDVSQVIRNMSSEFY  
ATQLRARITDETHPAAYEPEFYLQDDGGTAHLSAVSEDGSAVAATSTINLYF  
GSKVLSRVSGILFNDEMDDFSSPNFINQFRVAPSPANFIKPGKQPLSSMCPSII  
LDKDGQVRMVVGASGGTQITTSVALAIINSLWFGYDVKRAVEEPRLHNQLLPNT  
TTVEKDIDQVVTAGLKIRHHHTVTP

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:**N-terminal Tags, His-tag and T7-tag

**Subcellular Location:** Cell membrane.

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:**PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.9

**Predicted Molecular Mass:** 30.3kDa

**Accurate Molecular Mass:** 34kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

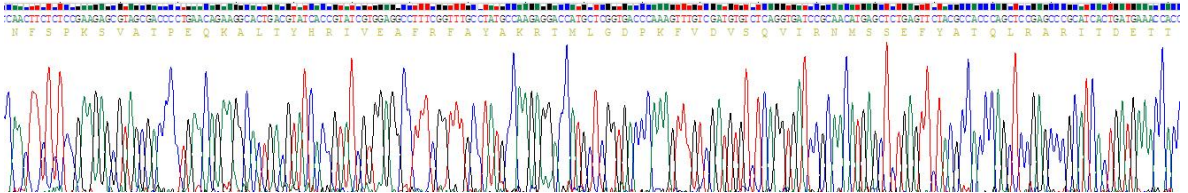


Figure 1. Gene Sequencing (Extract)

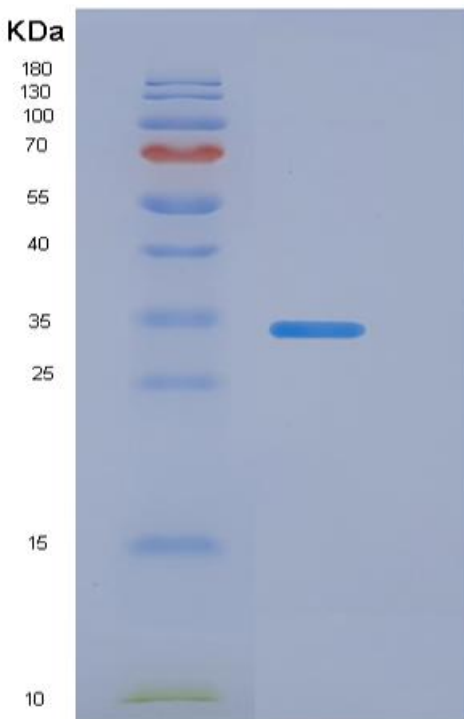


Figure 2. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any



issue if the kit was used in clinical diagnostic or any other procedures.