



Recombinant Tumor Necrosis Factor Receptor Superfamily, Member 8 (TNFRSF8)

Catalog No.: **TP04257** **50µg**

Sequence Information

Species: Human

Gene ID:943

Swiss Prot:P28908

Synonyms:CD30; KI-1; CD30L receptor; Ki-1 antigen; Lymphocyte activation antigen CD30

Residues:Met464~Lys595

MSQPLMETCHSVGAAYLESPLQDASPAGGPSSPRDLPEPRVSTEHTNKNKIEKI

YIMKADTVIVGTVKAELPEGRGLAGPAEPELEEELEADHTPHYPEQETEPPLGS

CSDVMLSVEEEGKEDPLPTAASGK

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags:N-terminal His Tag

Subcellular Location: Membrane, Cytoplasm.

Purity: >97%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 300µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.1

Predicted Molecular Mass: 17.8kDa

Accurate Molecular Mass: 21kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

CGAATTGGGATTCGGGACCAATTTGCTGTCCACAGTCATGCTAGCCATATGGCTGCCCGGGGGACAGGCGCTGCTGTGATGATGATGATGATGGCTGCTGCCATGGTATATCTCCTCTTAAAGTTAAACAAAATTATTCTAGAGGGGAAATGTTATCCGCTCACAAATTCGCCATATAGTGATCSTATTAA

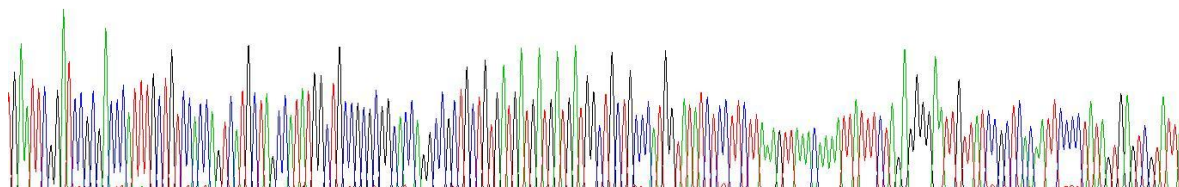


Figure 1. Gene Sequencing (Extract)

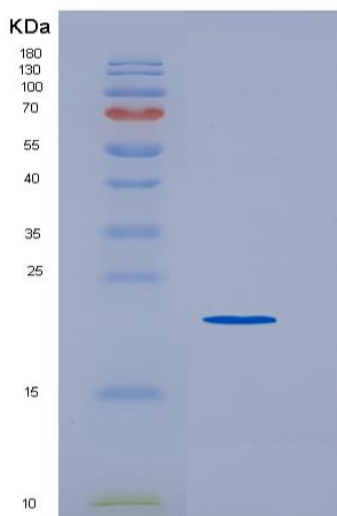


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.