



Recombinant Receptor Activator Of Nuclear Factor Kappa B Ligand (RANKL)

Catalog No.: TP04260 50µg

Sequence Information

Species: Human

Gene ID:8600

Swiss Prot:Q5T9Y4

Synonyms:CD254; TNFSF11; ODF; OPGL; TRANCE; HRANKL2; SOf; Tumor Necrosis Factor(ligand)superfamily Member 11; TNF-related activation-induced cytokine

Residues:Ile79~Ile247

ISEDGTHCIYRILRLHENADFQDTTLESQDTKLIPDSCRRIKQAFQGA VQKELQ
HIVGSQH IRAEKAMVDG SWLDLAKRSKLEAQPFAHLTINATDIPSGSHK VSLSS
WYHDRGWAKISNMTFSNGKLI VNOQDGFYLYANICFRHHETSGDLATEY LQLMV
YVTKTSI

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags:N-terminal His Tag

Subcellular Location:Membrane, Secreted, Cytoplasm, Chromosome.

Purity: >97%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 100µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.3

Predicted Molecular Mass: 22.9kDa

Accurate Molecular Mass: 22kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

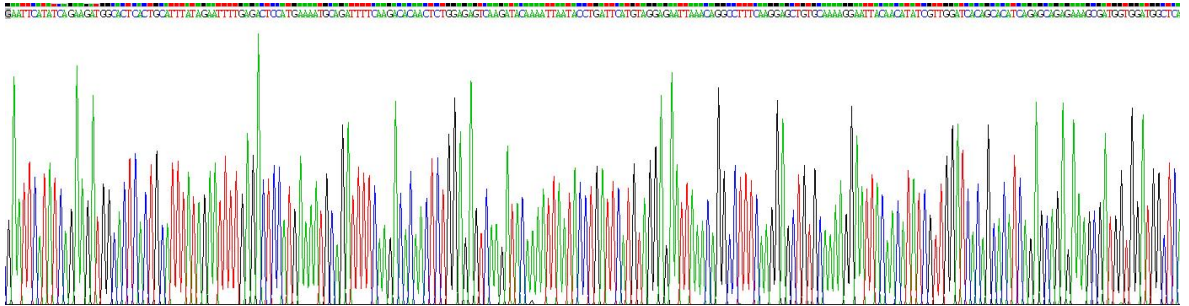


Figure 1. Gene Sequencing (Extract)

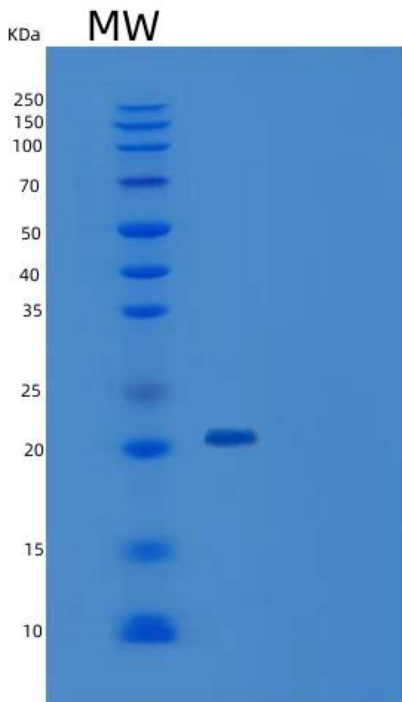


Figure 2. SDS-PAGE



[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.