

Recombinant Fibroblast Growth Factor 2, Basic (FGF2)

Catalog No.: TP04643

50µg

Sequence Information

Species: mouse

Gene ID:14173

Swiss Prot:P15655

Synonyms:B-FGF; BFGF; FGFB; HBGH-2;
Basic Fibroblast Growth Factor;
Heparin-binding growth factor 2

Residues:Pro10~Ser154

PALPEDGGAAFPPGHFKDKPKRLRYCKNGGFFLRIHPDGRVDGVREKSDPHVKLQL

QAEERGVSIVKGVCANRYLAMKEDGRLLASKCVTEECFFFERLESNNYNTYRSR

KYSSWYVALKRTGQYKLGSKTGPQKAILFLPMSAKS

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags:N-terminal His Tag

Subcellular Location: Nucleus, Secreted.

Purity: >97%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.7

Predicted Molecular Mass: 20.0kDa

Accurate Molecular Mass: 20kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was

determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

FALPEDGGAAFPFPCHFKDPKRLYCKNGGFFLRIHFDGRVDGVREKSDPHVKLQLQAEREGVYSTKGVCANRYLANKEDGRLLASKCVTEECFFFER

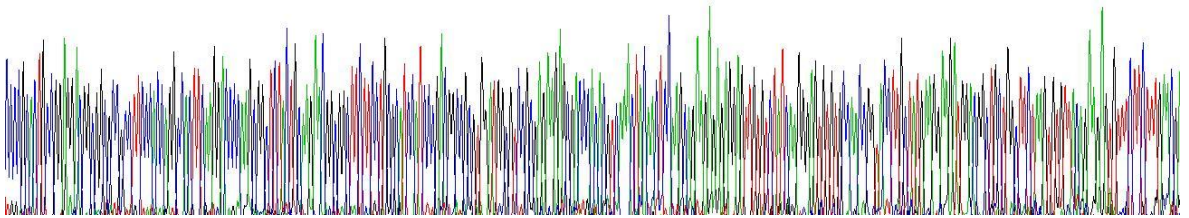


Figure 1. Gene Sequencing (Extract)

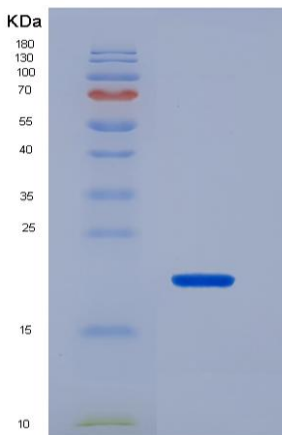


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.