

## Eukaryotic E-cadherin

Catalog No.: TP06402

50µg

### Sequence Information

**Species:** Human

**Gene ID:**999

**Swiss Prot:**P12830

**Synonyms:**CDHE; CD324; CDH1; Arc1; ECAD;  
CDH-E; E-CAD; LCAM; UVO;  
Cadherin 1 Type 1; Uvomorulin;  
E-Cadherin; Calcium-Dependent  
Adhesion Protein, Epithelial; Liver Cell  
Adhesion Molecule

**Residues:**Asp155~Ile707

DWVIPPISCPENEKGFPPKNLVQIKSNKDKGKVFYSITGQGADTPPVGVFIIE  
RETGWLVKTEPLDRERIATYTLFSAVSSNGNAVEDPMEILITVTDQNDNKPEF  
TQEVFKGSVMEGALPGTSVMEVTATDADDDVNTYNAAIAYTILSQDPELPDKNM  
FTINRNTGVISVTTGLDRESFPTYTLVVQAADLQGEGLSTTATAVITVTDTND  
NPPIFNPTTYKQVPENEANVVITTLKVTADAPNTPAWEAVYTI LNDDGGQFV  
VTTNPVNNDGILKTAKGLDFEAKQQYILHVAVTNVPFEVSLTTSTATVTVDVL  
DVNEAPIFVPPEKRVEVSEDFGVGQEITSYTAQEPDTFMEQKITYRIWRDTANW  
LEINPDTGAISTRAELDREDFEHVKNSTYTALIIATDNGSPVATGTGTL LLLILS  
DVNDNAPIPEPRTIFFCERNPKPQVINIIDADLPNTPSPFTAELTHGASANWTI  
QYNDPTQESIILKPKMALEVGDYKINLKLMDNQNKDQVTTLEVSVCDCEGAAGV  
CRKAQPVEAGLQI

### Product Information

**Source:** Eukaryotic expression.

**Host:** 293F cell

**Tags:** N-terminal His-Tag

**Subcellular Location:** Membrane, Golgi apparatus.

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 4.0

**Predicted Molecular Mass:** 62.0kDa

**Accurate Molecular Mass:** 75kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

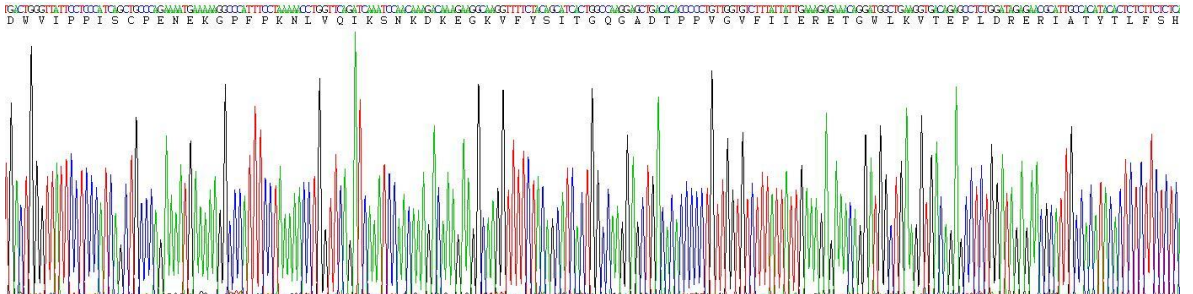


Figure 1. Gene Sequencing (Extract)

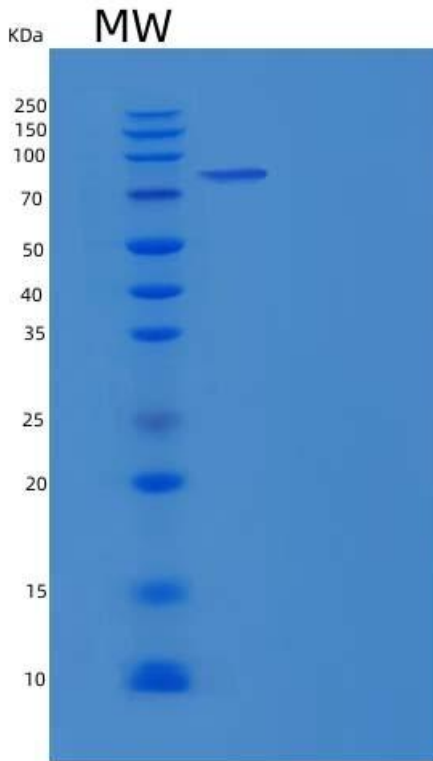


Figure 2. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.