

Eukaryotic Fc Fragment Of IgG Low Affinity IIIb Receptor (FcγR3B)

Catalog No.: TP05169 50μg

Sequence Information

Species: Human

Gene ID:2215

Swiss Prot:O75015

Synonyms:CD16b; CD16-B; FCGR3-B; FCG3;
FcγRIII; Low affinity immunoglobulin
gamma Fc region receptor III-B

Residues:Pro24~Ser203

PKAVVFLEPQWYSVLEKDSVTLKQCQGAYSPEDNSTQWFHNESLISSQASSYFID
AATVNDSGEYRCQTNLSTLSDPVQLEVHIGWLLLQAPRWVFKEDPIHLRCHSW
KNTALHKVTYLQNGKDRKYFHNSDFHHPKATLKDSGSYFCRGLVGSKNVSSET
VNITITQGLAVSTISSFS

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: N-terminal His-Tag

Subcellular Location: Membrane, Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 300μg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.9

Predicted Molecular Mass: 21.2kDa

Accurate Molecular Mass: 38-44kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

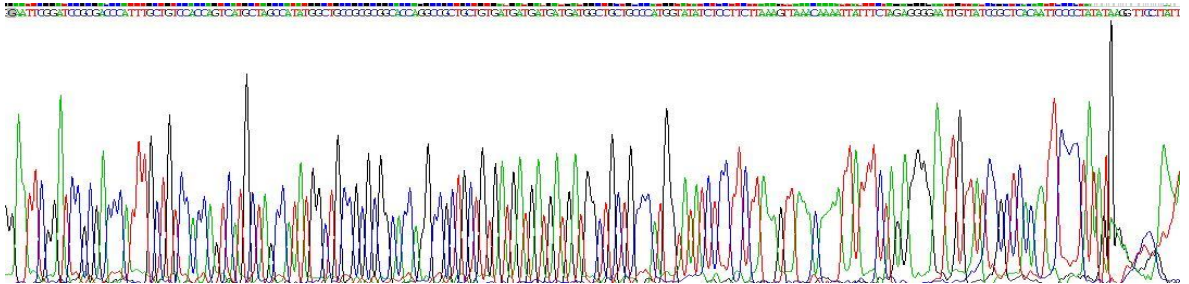


Figure 1. Gene Sequencing (Extract)

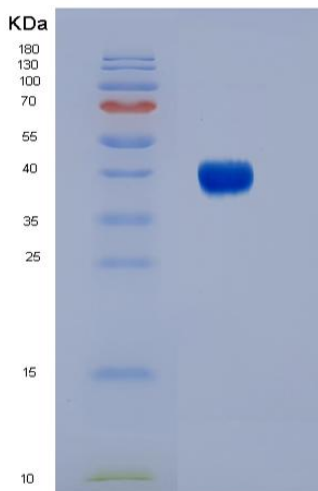


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.