



Recombinant Vascular Endothelial Growth Factor Receptor 1 (VEGFR1)

Catalog No.: TP04956

50µg

Sequence Information

Species: Human

Gene ID:2321

Swiss Prot:P17948

Synonyms:Flt1; Fms-Related Tyrosine Kinase 1;
Vascular Permeability Factor Receptor;
Fms-like tyrosine kinase 1;
Tyrosine-protein kinase receptor FLT

Residues:Ser27~Ile328

SKLKDP ELSLKG TQHIMQAGQTLHLQCRGEAAHKWSLP E M VSKESERLSITKSA
CGRNGKQFCSTLT LNTAQANHTGFYSCKYLAVPTSKKKETESA IYIFISDTGRP
FVEMYSEIPEIIHMTGRELVIPCRVTSPNITVTLKKFPLDTLIPDGKRIIWDS
RKGFIISNATYKEIGLLTCEATVNGHLYKTNYLTHRQTNTIIDVQISTPRPVKL
LRGHTLVLNCTATTPLNTRVQMTWSYPDEKNKRASVRRRIDQSN SHANIFYSVL
TIDKMQNKDKGLYTCRVRSGPSFKSVNTSVHI

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags:N-terminal His Tag

Subcellular Location:Membrane.

Purity: >97%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 350µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 10.0

Predicted Molecular Mass: 35.7kDa

Accurate Molecular Mass: 50kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

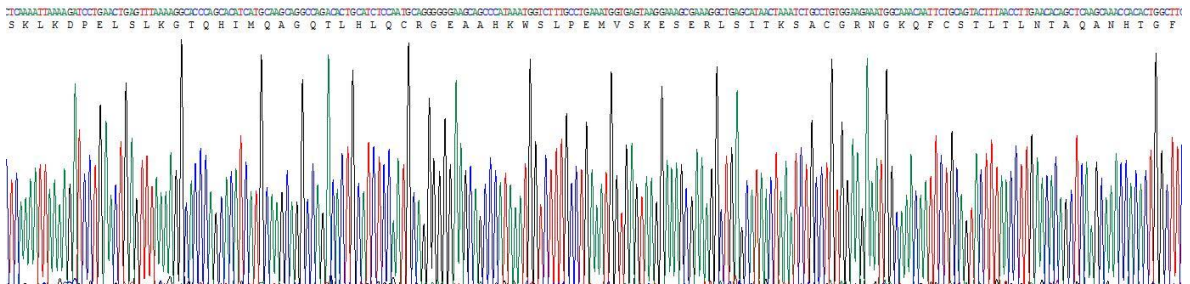


Figure 1. Gene Sequencing (Extract)

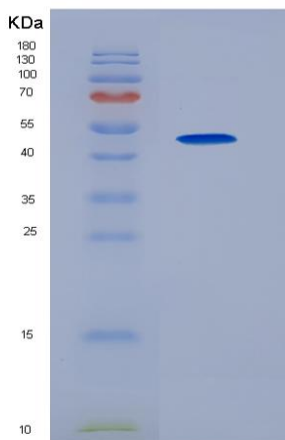


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.