



## Recombinant Growth Factor Receptor Bound Protein 2 (Grb2)

Catalog No.: TP05676

50µg

### Sequence Information

**Species:** Human

**Gene ID:**2885

**Swiss Prot:**P62993

**Synonyms:**ASH; NCKAP2; EGFRBP-GRB2;

Grb3-3; Adapter protein GRB2; Protein

Ash; SH2/SH3 adapter GRB2

**Residues:**Met 1-Val217

MEAIAKYDFKATADDELSFKRGDILKVLNEECDQNWYKAELNGKDGFIKPNYIE

MKPHPWFFGKIPRAKAEEMLSKQRHDGAF LIRESSESAPGDFSLSVKFGNDVQHF

KVLRDGAGKYFLWVVKFNSLNELVDYHRSTSVSRNQQIFLRDIEQVPQQPTYVQ

ALFDFDPQEDGELGFRRGDFIHVMDNSDPNWWKGACHGQTGMFPRNYVTPVNRN

V

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:**N-terminal His Tag

**Subcellular Location:**Nucleus; Cytoplasm;

Endosome; Golgi apparatus.

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:**PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.5

**Predicted Molecular Mass:** 29.1kDa

**Accurate Molecular Mass:** 28kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

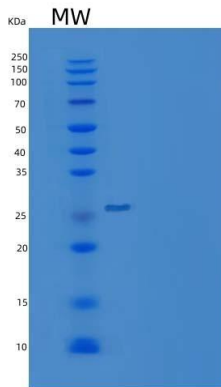


Figure 2. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.