



Recombinant Heat Shock 70kDa Protein 1B (HSPA1B)

Catalog No.: TP08686

10µg

Sequence Information

Species: Mouse

Gene ID:15511

Swiss Prot:P17879

Synonyms:Hsp70a1

Residues:Met1-Asp642

MAKNTAIGIDLGTTYSCVGVFQHGKVEIIANDQGNRTTPSYVAFTDTERLIGDAAKNQVALNPQNTVFDKRLIGRKFQDAVVQSDMKH
WPFQVVNDGDGPKVQVNYKGESRSFFPEEISSMVLTKMKEIAEAYLGHPVTNAVITVPAYFNDSQRQATKDAGVIAGLNLVRIINEPTA
AAIAYGLDRTGKGERNVLIIFDLGGGTFDVSILTIDDGIFEVKATAGDTHLGGEDFDNRLVSHFVEEFKRKHKKDISQNKRAVRRRLRTAC
ERAKRTLSSSTQASLEIDSLFEGIDFYTSITRARFEELCSDLFRGTLEPVEKALRDAKMDKAQIHDLVLVGGSTRIPKVQKLLQDFNNG
RDLNKSINPDEAVAYGAAVQAAILMGDKSENVQDLLLLDVAPLSLGLLETAGGVMTALIKRNSTIPTKQTQTFTTYSDNQPGVLIQVYEG
ERAMTRDNNLLGRFELS GIPAPRGVPQIEVTFDIDANGILNVTATDKSTGKANKITITNDKGRLSKEEIERMVQEAERYKAEDEVQRD
RVAAKNALESYAFNMKSAVEDEGLKGLSEADKKKVLDKCQEVISWLDSNTLADKEEFVHKREELERVCSPIISGLYQGAGAPGAGGFG
AQAPPKGASGSGPTIEEVD

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His-Tag.

Subcellular Location: Membrane, Cytoplasm, Golgi apparatus.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.1

Predicted Molecular Mass: 70.1kDa

Accurate Molecular Mass: 70kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

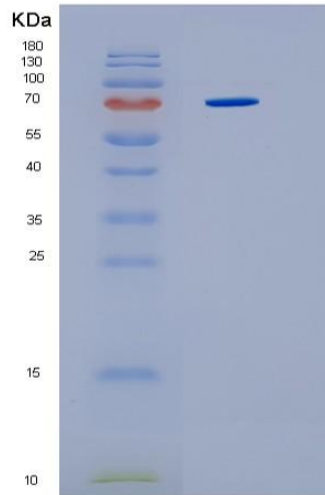


Figure . SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.