

Eukaryotic Interleukin 13 (IL13)

Catalog No.: TP05869

50µg

Sequence Information

Species: Human

Gene ID:3596

Swiss Prot:P35225

Synonyms:ALRH; BHR1; P600; IL-13; NC30

Residues:Leu25-Asn146

LTCLGGFASPGPVPPSTALRELIEELVNITQNQKAPLCNGSMVWSINLTAGMYC

AALESLINVSGCSAIEKTQRMMSGFCPHKVSAGQFSSLHVRDTKIEVAQFVKDL

LLHLKKLFREGRFN

Product Information

Source: Eukaryotic expression.

Host: *E. coli*

Tags:N-terminal His Tag and C-terminal Fc Region of Human IgG1

Subcellular Location:Secreted.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation:PBS, pH7.4, containing 0.1% SKL, 5% Trehalose.

Original Concentration: 80µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.2

Predicted Molecular Mass: 40.5kDa

Accurate Molecular Mass: 68kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is

less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

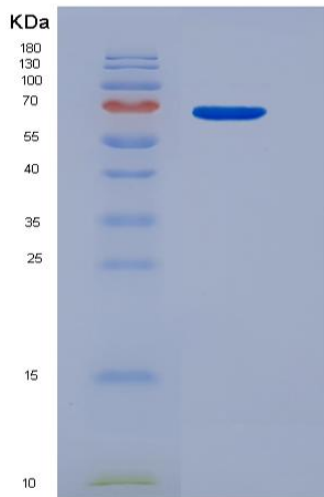


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.