



## Recombinant Human Plasminogen (Plg) Protein

Catalog No.: TP08983

50µg

### Sequence Information

**Species:** Human

**Gene ID:**5340

**Swiss Prot:**P00747

**Synonyms:**PL; Plasmin

**Residues:**Glu20-Asn810

EPLDDYVNTQGASLFSVTKKQLGAGSIEECAAKCEEDEEFTCRAFQYHSKEQQCVIMAENRKSSIIIRM  
RDVVLFEKKVYLSECKTGNKKNYRGTMSKTKNGITCQKWSSTSPHRPRFSPATHPSEGLEENYCRNPDN  
DPQGPWCYTTDPEKRYDYCDILECEEECMHCSGENYDGKISKTMGLECQAWDSQSPHAHGYIPSKFPN  
KNLKKNYCRNPDRELRPWCFTTDPNKRWELCDIPRCTTPPPSSGPTYQCLKGTGENYRGNVAVTVSGHT  
CQHWSAQTPHTHNRTPENFPCKNLDENYCRNPDGKRAPWCHTTNSQVRWEYCKIPSCDSSPVSTEQLAP  
TAPPELTPEVQDCYHGDGQSYRGTSSTTTTGKKCQSWSSMTPHRHQKTPENYPNAGLTMNYCRNPDADK  
GPWCFTTDPVSRWEYCNLKKCSGTEASVVAPPPVLLPDVETPSEEDCMFNGKGYRGRATTVTGTGTPC  
QDWAAQEPHRHSIFTPETNPRAGLEKNYCRNPDGVDGGPWCYTTNPRKLYDYCDVPQCAAPSFDCGKPKQ  
VEPKKCPGRVVGGCVAPHSWPWQVSLRTRFGMHFCGGTLISPEWVLTAAHCLEKSPRPSSYKVIILGAH  
QEVNLEPHVQEIIEVSRLFLEPTRKDIALLLKLSPPAVITDKVIPACLSPNYVADRTECFITGWGETQG  
TFGAGLLKEAQLPVIENKVCNRYEFLNGRVQSTELCAGHLAGGTDSCQGDSGGPLVCFEKDKYIILQGV  
SWGLGCARPKNKPGVYVRSRFTWIEGVMRNN

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Membrane.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS PH=7.2, 5% Trehalose.

**Original Concentration:** 400µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 8.5

**Predicted Molecular Mass:** 91.7kDa

**Accurate Molecular Mass:** 92kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in PBS (pH7.2) or ddH<sub>2</sub>O to a concentration of 0.1-0.5mg/mL. Do not vortex.

### [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at  $-80^{\circ}\text{C}$  for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at  $37^{\circ}\text{C}$  for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

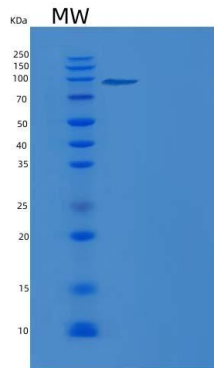


Figure . SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.