

Recombinant Cluster Of Differentiation 55 (CD55)

Catalog No.: TP08748

50µg

Sequence Information

Species: Rat

Gene ID:

Swiss Prot:A0A8I5XWT0

Synonyms:DAF; CR; TC; Complement
decay-accelerating factor; Decay
Accelerating Factor

Residues:Pro254~Glu372

DQGFILVGNSSVHCTLNGEVGEWSSAPPKCIERHKVPTKTPPEVDIPSTMPHES
TLINFPSTRAPLSHKPTTVKVPGRVKPTSHKPTEVKVPGTRVPPTSQKHTTVN
VPSTRVPATSQ

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 10.3

Predicted Molecular Mass: 14.7kDa

Accurate Molecular Mass: 16kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

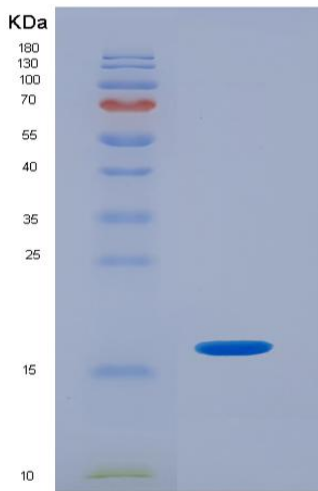


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.