

Recombinant Semaphorin 4D (SEMA4D)

Catalog No.: TP08750

50µg

Sequence Information

Species: Rat

Gene ID:

Swiss Prot:A0A8I6GJZ1

Synonyms:CD100; A8; BB18; GR3; M-Sema G;

M-Sema-G; SEMAJ; Coll-4; Sema

Domain,Immunoglobulin

Domain(Ig),Transmembrane

Domain(TM)and Short Cytoplasmic

Domain 4D

Residues:Leu145~Gly494

LGKSEDKGRCRPFDPAHSYTSVMVGGELYSGTSYNFLGSEPIISRNSSHSPLRT
EYAIPLWNEPSFVFADVIHKSPDGTEAEDDKVYFFFTVEVSVEYEFVFKLMIPRV
ARVCKGDQGGRLTLQKKWTSFLKARLICSRPDSGLVFNILQDVFVLRAPGLKEP
VFYAVFTPQLNNGVLSAVCAAYTLSTVEAVFSRGKYMQSATVEQSHTKWVRYNGP
VPTPRPGACIDSEARAANYTSSLNLPDKTLQFVKDHPLMDDSVTPIDNRPKLIK
KDVNYTQIVVDRITQALDGTFFYDVMFLSTDRGALHKAVILAKEVHVVEETQLFQD
FEPVLSLLLSSKKGRKFVYAGSNSG

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and GST Tag

Subcellular Location: Cytoplasm, Extracellular matrix.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.2

Predicted Molecular Mass: 68.9kDa

Accurate Molecular Mass: 69kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

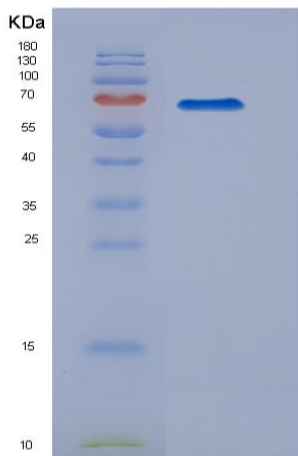


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.