



Recombinant Human Microtubule Associated Protein 1 Light Chain 3 Beta 2 (MAP1LC3b2)

Catalog No.: TP02895

50µg

Sequence Information

Species: Human

Gene ID:643246

Swiss Prot:A6NCE7

Synonyms:Microtubule-associated proteins

1A/1B light chain 3 beta 2, ATG8G

Residues:Met1~Gly120

MPSEKTFKQRRTFEQRVEDVRLIREQHPTKIPVIIERYKGKQLPVLDKTKFLV

PDHVNMSSELIKIIRRLQLNANQAFFLLVNGHSMVSVSTPISEVYESEKDEDEGF

LYMVCASQETFG

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.5

Predicted Molecular Mass: 17.7kDa

Accurate Molecular Mass: 20kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

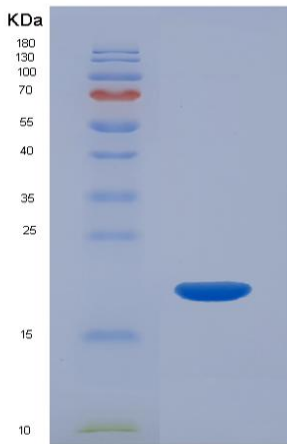


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.