

Recombinant Serine/threonine-protein kinase mTOR (mTOR)

Catalog No.: TP08789

50µg

Sequence Information

Species: Human

Gene ID:2475

Swiss Prot:P42345

Synonyms:FRAP; FRAP1; FRAP2; MTOR;

RAFT1; RAPT1; FK506 Binding Protein

12 Rapamycin Associated Protein;

Mammalian Target Of Rapamycin;

FKBP12-Rapamycin

Complex-Associated Protein 1;

Rapamycin and FKBP12 target

Residues:Ala2226~Val2488

AVIPLSTNSGLIGWVPHCDTLHALIRDYREKKKILLNIEHRIMLRMAPDYDHLT
LMQKVEVFEHAVNNTAGDDLAKLLWLKSPSEVWFDRRTNYTRSLAVMSMVGVI
LGLGDRHPSNMLDRLSGKILHIDFGDCFEVAMTREKFPEKIPFRLTRMLTNAM
EVTGLDGNYRITCHTVMEVLREHKDSVMAVLEAFVYDPLLNWRLMDTNTKGNKR
SRTRDSYSAGQSVEILDGVELGEPAHKKTGTTVPESIHSEFIGDGLV

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Nucleus, Mitochondrion, Lysosome, Cytoplasm, Golgi apparatus, Endoplasmic reticulum lumen.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 500µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.1

Predicted Molecular Mass: 33.7kDa

Accurate Molecular Mass: 36kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

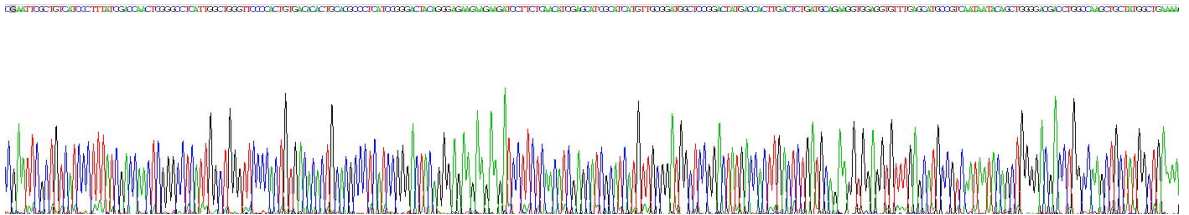


Figure 1. Gene Sequencing (Extract)

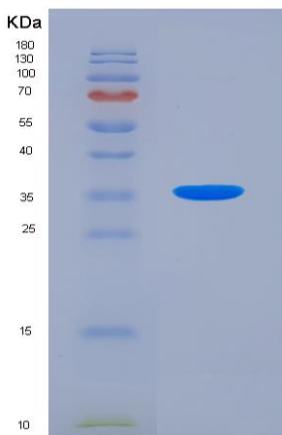


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

