

Recombinant Receptor Activator Of Nuclear Factor Kappa B Ligand (RANKL)

Catalog No.: TP06065 50µg

Sequence Information

Species: Human

Gene ID:8600

Swiss Prot:O14788

Synonyms:CD254; ODF; OPGL; RANKL;

TNFSF11; CD254; Osteoclast

differentiation factor; Receptor activator

of nuclear factor kappa-B ligand; tumor

necrosis factor ligand superfamily

member 11

Residues:Ser63~Ile247

SVALFFYFRAQMDPNRISEGDGTHCIYRILRLHENADFQDTTLESQDTKLIPDSC

RRIKQAFQGAQVQKELQHIVGSQHIRAEKAMVDGSLDLAKRSKLEAQPFAHLTI

NATDIPSGSHKVSLSWYHDRGWAKISNMTFSNGKLIVNQDGFYLYANICFRH

HETSGDLATEYLQLMVYVTKTSI

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and GST Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 500µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.7

Predicted Molecular Mass: 51.0kDa

Accurate Molecular Mass: 51kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

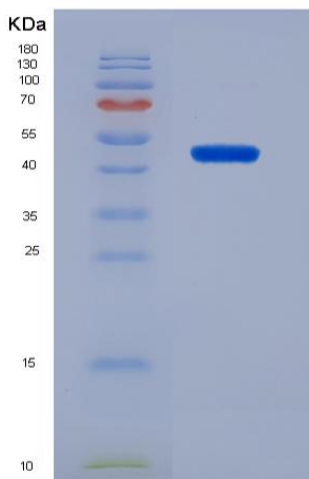


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.