

## Recombinant Human Ectonucleoside triphosphate diphosphohydrolase 3 (ENTPD3)

Catalog No.: TP02022

50µg

### Sequence Information

**Species:** Human

**Gene ID:**956

**Swiss Prot:**O75355

**Synonyms:**Ectonucleoside triphosphate  
diphosphohydrolase 3, CD39L3, HB6,  
NTPDase-3

**Residues:**Gln44-Pro485

QIHKQEVLPPLGKYGIVLDAGSSRTTVYVYQWPAEKENNTGVVSQTFKCSVKGS  
GISSYGNNPQDVPRAFECEMQKVKGQVPSHLHGSTPIHLGATAGMRLRLQNET  
AANEVLESIQSYFKSQPFDFRGAQIIISGQEEGVYGWITANYLMGNFLEKNLWHM  
WVHPHGVETTALDLGGASTQISFVAGEKMDLNTSDIMQVSLYGYVYTLYTHSF  
QCYGRNEAEKKFLAMLLQNSPTKNHLTNP CYPRDYSISFTMGHVFDSLCTVDQR  
PESYNPNVDVITFEGTGDPSLCCKEKVASIFDFKACHDQETCSFDGVYQPKIKGPF  
VAFAGFYTASALNLSGSFSLDTFNSSTWNFCSQNWSQLPLLLPKFDEVYARSY  
CFSANYIYHLFVNGYKFTEETWPQIHFEKEVGNSSIAWSLGMYLSLTNQIPAES  
PLIRLPIEPP

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Cell membrane

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:** 52.0kDa

**Accurate Molecular Mass:** 52kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

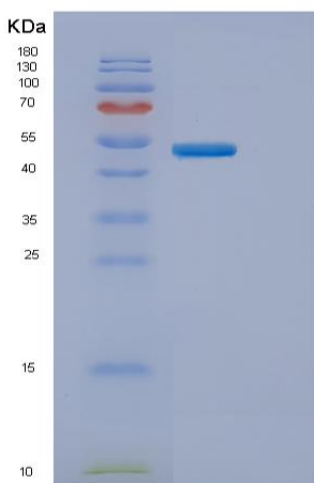


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.