

Recombinant Human Cartilage-associated protein (CRTAP)

Catalog No.: TP01721

50µg

Sequence Information

Species: Human

Gene ID:10491

Swiss Prot:O75718

Synonyms:Cartilage-associated protein
precursor , CASP, LEPREL3, OI7

Residues:Gln27-Ser401

QYERYSFERSFPRDELMPLESAYRHALDKYSGEHWAESVGYLEISLRHLRLLRDS
EAFCHRNCSAAPQPEPAAGLASYPELRLFGGLLRRAHCLKRCKQGLPAFRQSQP
SREVLADFQRREPYKFLQFAYFKANNLPKAIAAAHTFLLKHPDDEMKNRMAYY
KSLPGAEDYIKDLETKSYESLFIKAVRAYNGENWRTSITDMELALPDFFKAFYE
CLAACEGSREIKDFKDFYLSIADHYVEVLECKIQCEENLTPVIGGYPVEKFVAT
MYHYLQFAYYKLNLDKNAAPCAVSYLLFDQNDKVMQQNLVYYQYHRDTWGLSDE
HFQPRPEAVQFFNVTTLQKELYDFAKENIMDDDEGEVVEYVDDLLELEETS

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted, extracellular space, extracellular matrix.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.1

Predicted Molecular Mass: 46.4kDa

Accurate Molecular Mass: 46kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

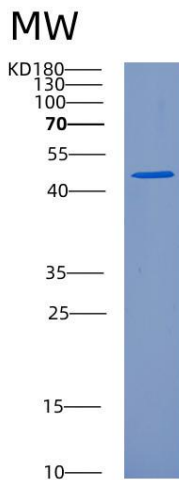


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.