



Recombinant Human Growth Arrest And DNA Damage Inducible Protein Gamma (GADD45g)

Catalog No.: **TP06135** **50µg**

Sequence Information

Species: Human

Gene ID:10912

Swiss Prot:O95257

Synonyms:Growth Arrest and DNA

Damage-Inducible Protein GADD45

Gamma; Cytokine-Responsive Protein

CR6; DNA Damage-Inducible Transcript

2 Protein; DDIT-2; GADD45G; CR6;

DDIT2

Residues:Met1-Glu159

MTLEEVRGQDTVPESTARMQGAGKALHELLLSAQRQGCLTAGVYESAKVLNVDP

DNVTFCVLAAGEEDEGDIALQIHFTLIQAFCCENDIDIVRVGDVQRLAAIVGAG

EEAGAPGDLHCILISNPNEAWKDPALKLSLFCEESRSVNDWVPSITLPE

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Nucleus, Cytoplasm

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.0

Predicted Molecular Mass: 21.0kDa

Accurate Molecular Mass: 24kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

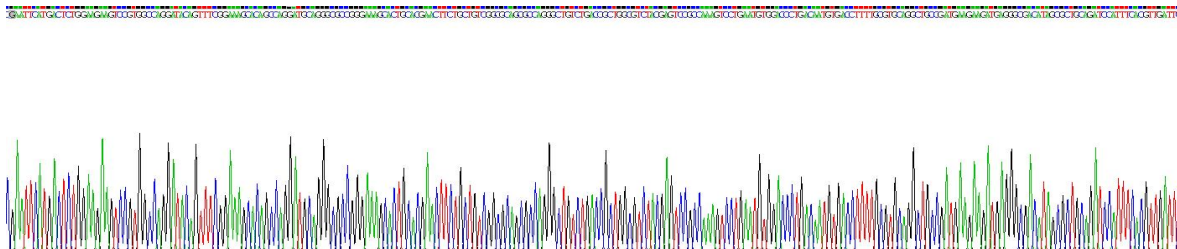


Figure 1. Gene Sequencing (Extract)

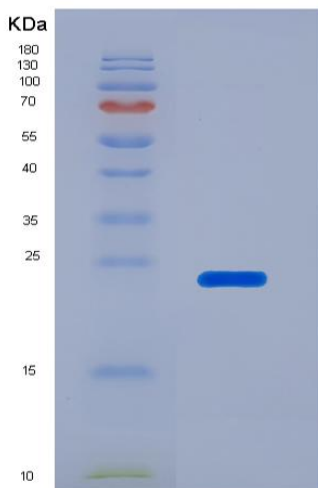


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

