

## Recombinant Human Antithrombin (AT)

Catalog No.: TP03902

50µg

### Sequence Information

**Species:** Human

**Gene ID:**462

**Swiss Prot:**P01008

**Synonyms:**AT3; ATIII; SERPINC1;

Anti-Thrombin Antibodies; Serpin

Peptidase Inhibitor Clade C Member 1;

Coding Sequence Signal Peptide

Antithrombin Part 1

**Residues:**His33~Lys464

HGSPVDICTAKPRDIPMNPNCIYRSPEKKATEDEGSEQKIPEATNRRVWELSKA  
NSRFATTFYQHLADSKNDNDNIFLSPLSISTAFAMTKLGACNDTLQQLMEVFKF  
DTISEKTSQIHHFFAKLNCRLYRKANKSSKLVSANRLF GDKSLTFNETYQDIS  
ELVYGAKLQPLDFKENAEQSRAAINKWVSNKTEGRITDVIPSEAINELTVLVLV  
NTIYFKGLWKS KFS PENTRKELFYKADGESCSASMMYQEGKFRYRRVAEGTQVL  
ELPFKGDDITMVLILPKPEKSLAKVEKELTPEVLQEWLDELEEMMLVVHMPFRF  
IEDGFS LKEQLQDMGLVDLFSPEKSKLPGIVAEGRDDLYVSDAFHKAFLEVNEE  
GSEAAASTAVVIAGRSLNPNRVTFKANRPFLVFIREVPLNTIIFMGRVANPCVK

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His Tag

**Subcellular Location:** Secreted, Extracellular matrix

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.2

**Predicted Molecular Mass:** 50.3kDa

**Accurate Molecular Mass:** 50kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

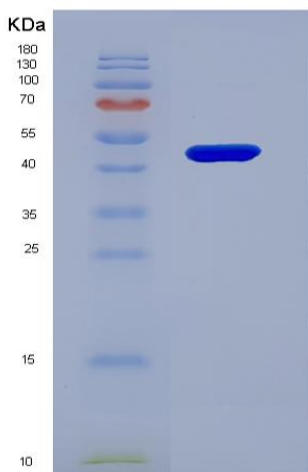


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.