

Recombinant Human Alpha-1-Antitrypsin (a1AT)

Catalog No.: TP01175

50µg

Sequence Information

Species: Human

Gene ID:5265

Swiss Prot:P01009

Synonyms:Alpha-1 antiproteinase (Serine proteinase inhibitor), PI, A1A, AAT, PI1, A1AT, MGC9222, PRO2275, MGC23330, SERP1NA1.

Residues:Glu25~Lys418

EDPQGDAAQKTDTSHHDDQDHTFNKITPNLAFAFSLYRQLAHQSNSTNIFFSP
VSIATAFAMLSLGTKADTHDEILEGLNFNLTEIPEAQIHEGFQELLRTLNPDS
QLQLTTGNGLFLESEGLKLVDFLEDVKKLYHSEAFVNFVDTEEAKKQINDYVE
KGTQGKIVDLVKELDRDVFALVNYIFFKGKWERPFVVKDTEEDFHVDQVTTV
KVPMMKRLGMFNIQHCKKLSSWVLLMKYLGNATAIFFLPDEGKLQHLENELTHD
IITKFLNEDRRSASLHLPKLSITGTYDLKSVLGQLGITKVFSGADLSGVTEE
APLKLSKAVHKAVLTIDEKGTEAAGAMFLEAIPMSIPPEVKFNKPFVFLMIEQN
TKSPLFMGKVVNPTQK

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted, Endoplasmic reticulum lumen

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.4

Predicted Molecular Mass: 48.0kDa

Accurate Molecular Mass: 48kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

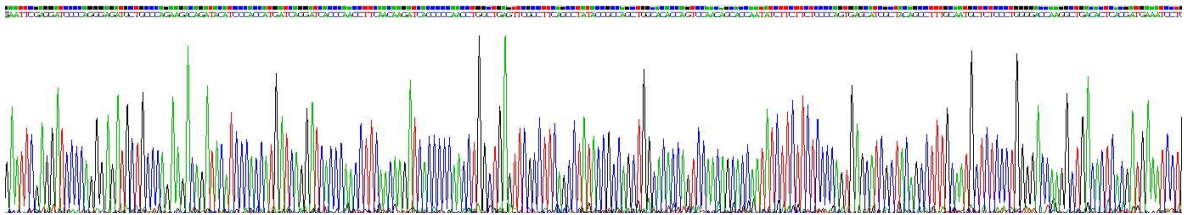


Figure 1. Gene Sequencing (Extract)

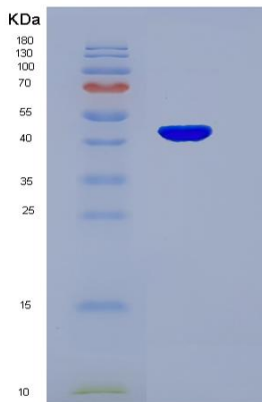


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.