



Recombinant Glycophorin A (GYPA)

Catalog No.: **TP09182** **100µg**

Sequence Information

Species: Human

Gene ID:2993

Swiss Prot:A0A2R8Y7F9

Synonyms:CD235a; MN; MNS;
GPA; GPErik; GPSAT; PAS-2;
GpMiIII; HGpMiIII; HGpMiV;
HGpMiX; HGpMiXI; HGpSta;
Sialoglycoprotein alpha; MN
sialoglycoprotein

Residues:Leu20~Glu91

LSTTEVAMHTSTSSSVTKSYISSQTNDTHKRDTYAATPRAHEVSEISVRV

YPPEETGER VQLAHHFSEP E

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His and GST Tag

Subcellular Location: Membrane, Chromosome

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS (PH7.4) , containing 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.3

Predicted Molecular Mass: 11.9kDa

Accurate Molecular Mass: 15kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage



ITAGTCTCGGGGAGTCCCTAACTCTCTGGGAGGTATCTATCGGAGGATCGGAGGTGGATATGCTCCCTAGCGCTGAGTGGAGGTTGTAGAGTTATCTGGAGGTGGGAGGTGTAATCTCATTCGAGAGG
LSTTEVAMHTSTSSSVTKSVISSQNDTHKRDTVAAATPRAHEVSEISVVRTVVPPEEETGERVOLAHHESEPP

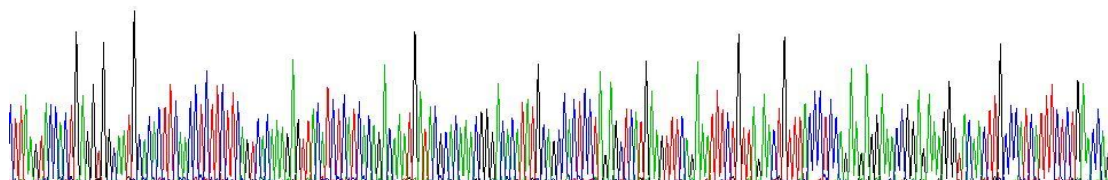


Figure 1. Gene Sequencing (Extract)

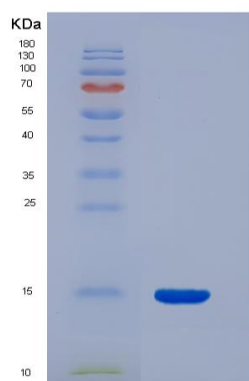


Figure 2. SDS-PAGE

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.