



## Recombinant Mouse Low Density Lipoprotein Receptor

### Related Protein 1 (LRP1)

Catalog No.: TP09271 1mg

#### Sequence Information

**Species:** Mouse

**Gene ID:**16971

**Swiss Prot:**Q91ZX7

**Synonyms:**CD91; LRP; A2MR; APOER;

APR; TGFBR5;

Alpha-2-Macroglobulin

Receptor; Prolow-density

lipoprotein receptor-related

protein 1; Apolipoprotein E

receptor

**Residues:**Ile692~Cys962

IFVTSKTVLWPNGLSLDIPAGRLYWDAFYDRIETILLNGTDRKIVYEGPELNHAFGLCHHGNYLFWTEYRSGSV  
YRLERGVAGAPPTVTLLRSERPPIFEIRMYDAQQQVGTNKRNVNNGGSSSLCLATPGSRQCACAEDQVLDTDGV  
TCLANPSYVPPPPQCQPGEFACANNRCIQERWKCDGDNDCLDNSDEAPALCHQHTCPSDRFKCENNRCIPNRWLCD  
GDNDCGNSEDESATCSARTCPPNQFSCASGRICIPISWTCDLDDDC

#### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** N-terminal His Tag

**Subcellular Location:** Chromosome.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS pH7.2, 0.01% sarcosyl, 5%Trehalose .

**Original Concentration:** 1mg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.4

**Predicted Molecular Mass:** 33.8kDa

**Accurate Molecular Mass:** 34kDa as determined by SDS-PAGE reducing conditions.

#### **[ USAGE ]**

Reconstitute in PBS pH7.2 or ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

#### **[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.




Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## **[ IDENTIFICATION ]**

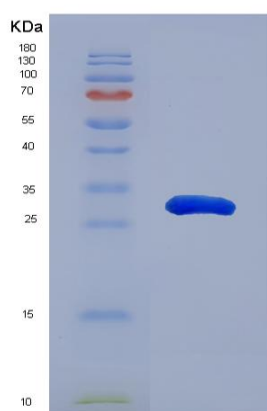

**Figure 1. SDS-PAGE**

## **[ IMPORTANT NOTE ]**

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.