

Recombinant Human Ras Homolog Gene Family, Member Q (RHOQ)

Catalog No.: TP03738

100µg

Sequence Information

Species: Human

Gene ID:23433

Swiss Prot:P17081

Synonyms:Rho-related GTP-binding protein

RhoQ precursor , ARHQ, RASL7A,

TC10, TC10A

Residues:Met1~Cys202

MAHGPGALMLKCVVVG DGAVGKTCLLMSYANDAFPEEVPTVFDHYAVSVTVGG

KQYLLGLYDTAGQEDYDRLRPLSYPM TDVFLICFSV VNPASFQNVKEEWPELK

EYAPNVPFLLIGTQIDL RDDPKTLARLNDMKEK PICVEQGQKLAK EIGACCYVE

CSALTQKGLKTVFDEAIIAILTPKHTVKKRIGSRCINCC

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His and GST Tag

Subcellular Location: Membrane, Cytoplasm

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.9

Predicted Molecular Mass: 52.3kDa

Accurate Molecular Mass: 52kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

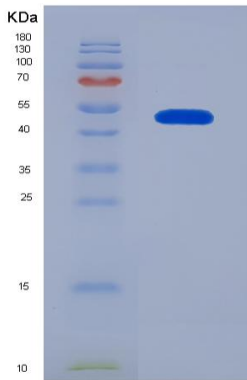


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.