



Polyclonal Antibody to Gonadotropin Releasing Hormone (GnRH)

PA12407 100ul

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Recombinant Protein affinity chromatography.

Traits: Liquid

Concentration: 500ug/ml

Ig Isotype:IgG

UOM: 100µl

Cross Reactivity: Human; Mouse .

Applications: IHC,ELISA.

[IMMUNOGEN]

Immunogen: Recombinant Human GnRH protein expressed in E.coli.

[APPLICATIONS]

IHC 1:100-1:300

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in 0.01M PBS, pH7.4, containing 0.05% Proclin-300, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 4°C for frequent use.

Aliquot and store at -20°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was

determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

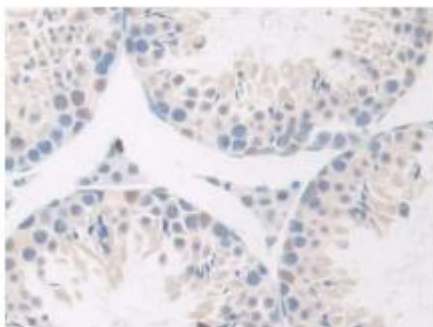

DAB staining on IHC-P;

Samples: Mouse Testis Tissue;

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.