



Recombinant Matrix Metalloproteinase 9 (MMP9)

Catalog No.: TP09331

100µg

Sequence Information

Species: Human

Gene ID:4318

Swiss Prot:P14780

Synonyms:CLG4B

Residues:Ala20-Asp707

APRQRQSTLVLPFGDLRTNLTDRQLAEEYLYRYGYTRVAEMRGESKSLGPALLLLQKQLSLPETGELDSATLKAMRTPR
CGVPDLGRFQTFEGDLKWHHHNITYWIQNYSEDLPRVIDDAFARAFALWSAVTPLTFTRVYSRDADIVIQFGVAEHGD
GYFPDGDGDL LAHAFPPGPGIQGDAHFDDELWSLKGKVVVPTFRGNADGAACHFPFIFEGRSYSACTTDGRSDGLPWC
STTANYDTDDRFGFCPSERLYTQDGNADGKPCQFPFIFQGGSYSACTTDGRSDGYRWCATTANYDRDKLFGFCPTRADS
TVMGGNSAGELCVFPFTFLGKEYSTCTSEGRGDGRLWCATTSNFSDSKKWGFCPDQGYSLFLVAAHEFGHALGLDHSSV
PEALMYPMYRFTEGPPLHKDDVNGIRHLYGPRPEPEPRPPTTTTPQPTAPPTVCPTGPPTVHPSERPTAGPTGPPSAGP
TGPPTAGPSTATTVPLSPVDDACNVNIFDAIAEIGNQLYLFDKGKYWRFSEGRGSRPQGPFLIADKWPALPRKLDVFE
ERLSKKLFFFSGRQVWVYTGASVLGPRRLDKLGLGADVAQVTGALRSGRGKMLLFSGRRLWRFVDVKAQMVDPRSASEVD
RMFPGVPLDTHDVFQYREKAYFCQDRFYWRVSSRSELNQVDQVG YVTYDILQCPED

Product Information

Source: Eukaryotic expression.

Host: 293F

Tags: N-terminal His-Tag.

Tissue Specificity:

Subcellular Location: Secreted, extracellular space, extracellular matrix.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, 5%Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.6

Predicted Molecular Mass: 80kDa

Accurate Molecular Mass: 95kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-0.5mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

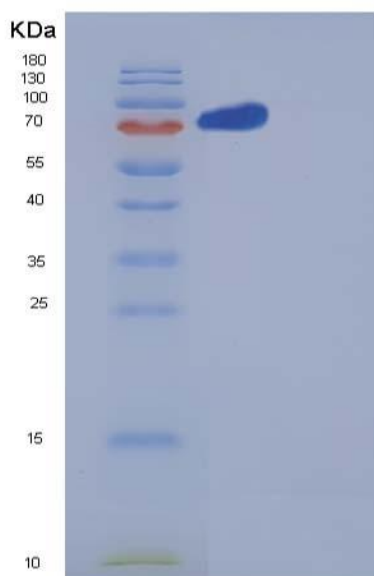


Figure 1. Recombinant Human MMP9 SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.