

Eukaryotic Calreticulin (CALR)

Catalog No.: TP07113

100µg

Sequence Information

Species: Human

Gene ID:811

Swiss Prot:P27797

Synonyms:cC1qR, CRT, FLJ26680, RO, SSA.

Residues:Glu18~Leu417

EPAVYFKEQFLDGDGWSRWIESKHKSDFGKFVLSSGKFYGDEEKDKGLQTSQD
ARFYALSASFEPFSNKGQTLVVQFTVKHEQNIDCGGGYVKLFPNSLDQTMHGD
SEYNIMFGPDICGPGTKKVHVIFNYKGKNVLINKDIRCKDDEFTHLYTLIVRPD
NTYEVKIDNSQVESGSLEDDWDFLPPKKIKDPDASKPEDWDERAKIDDPTDSKP
EDWDKPEHIPDPDAKKPEDWDEEMDGEWEPPIQNPEYKGEWKPRQIDNPDYKG
TWIHPEIDNPEYSPDPSIYAYDNFVGLGLDLWQVKSGTIFDNFLITNDEAYAE
FGNETWGVTKAAEKQMKDKQDEEQRLKEEEEDKKRKEEEEAEDEKEDDEDKDEDE
EDEEDKEEDEEDVPGQAKDEL

Product Information

Source: Eukaryotic expression.

Host: HEK293 Cells

Tags: N-terminal His-Tag

Subcellular Location: Secreted, Cytoplasm, Endoplasmic reticulum lumen, Extracellular matrix.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.0

Predicted Molecular Mass: 48.1kDa

Accurate Molecular Mass: 60kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

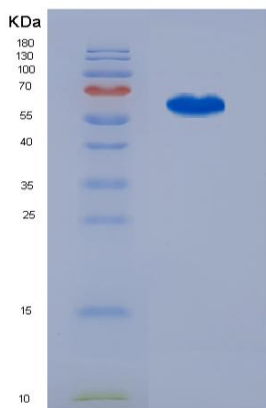


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.