



Recombinant Aldehyde Dehydrogenase 1 Family, Member A1 (ALDH1A1)

Catalog No.: TP09355

100µg

Sequence Information

Species: Mouse

Gene ID:11668

Swiss Prot:P24549

Synonyms:ALDH11; ALDH1; ALDC; PUMB1;

RALDH1; Retinaldehyde

Dehydrogenase 1; Aldehyde

dehydrogenase, cytosolic; Retinal

dehydrogenase 1

Residues:Met1~Arg326

MSSPAQPAVPAPLADLKIQHTKIFINNEWHNSVSGKKFPVLNPATEEVICHVEE

GDKADVDKAVKAARQAFQIGSPWRTMDASERGRLLNKLADLMERDRLLLLATMEA

LNCGKVFANAYLSDLGGCIKALKYCAGWADKIHGQTIPSDGDIITYTRREPIGV

CGQIIPWNFPMLMFIWKIGPALSCGNTVVVKPAEQTPLTALHLASLIKEAGFPP

GVVNIVPGYGPTAGAAISSHMDVDKVAFTGSTQVGKLIKEAAGKSNLKRVTLLEL

GGKSPCIVFADADLDIAVEFAHHGVFYHQGCCVAASRIFVEESVYDEFVKRSV

ER

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His Tag

Subcellular Location: Cytoplasm.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.75

Predicted Molecular Mass: 39.0kDa

Accurate Molecular Mass: 39kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

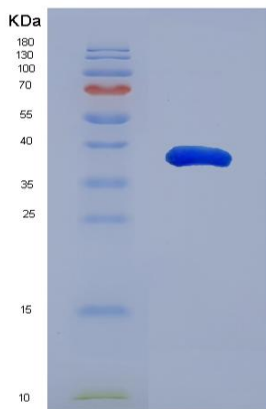


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.