

Recombinant Heat Shock Protein 75kDa, Mitochondrial (HSP75)

Catalog No.: TP09357

100µg

Sequence Information

Species: Human

Gene ID:10131

Swiss Prot:Q12931

Synonyms:TRAP1; HSP90L; TNF

Receptor-Associated Protein 1; Tumor
necrosis factor type 1 receptor-associated
protein

Residues:Ser60~His704

STQTAEDEKEPLHSIISSTESVQGSTSKHEFQAETKKLLDIVARSLYSEKEVFI
RELISNASDALEKLRHKLVS DGQALPEMEIHLQTNAEKGTTITIQDTGIGMTQEE
LVSNLGTIARSGSKAFLDALQNQAEASSKIIGQFGVGFYSAFMVADRVEVYSRS
AAPGSLGYQWLSDGSGVFEIAEASGVRTGTKIIHLKSDCKEFSSEARVRDVVT
KYSNFVSFPLYLNGRRMNTLQAIWMDPKDVREWQHEEFYRYVAQAHDKPRYTL
HYKTDAPLNIRSFYVPDMKPSMFDVSREL GSSVALYSRKVLIQTKATDILPKW
LRFIRGVVDSEDIPLNLSRELLQESALIRKLRDVLQQRLIKFFIDQSKKDAEKY
AKFFEDYGLFMREGIVTATEQEVKEDIAKLLRYESSALPSGQLTSLSEYASRMR
AGTRNIYYLCAPNRHLAEHSPYYEAMKKKDTEVLFCFEQFDELTLHLREFDKK
KLISVETDIVVDHYKEEFEDRSPAAECLSEKETEELMAWMNRNVLGSRVTNVKV
TLRLDTHPAMVTVLEMGAAHF LRMQQLAKTQEERAQLLQPTLEINPRHALIKK
LNQLRASEPGLAQLLVLDQIYENAMIAAGLVDDPRAMVGRLNELLVKALERH

Product Information

Source: Recombinant expression.

Host: *E.coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.5

Predicted Molecular Mass: 77.3kDa

Accurate Molecular Mass: 84kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

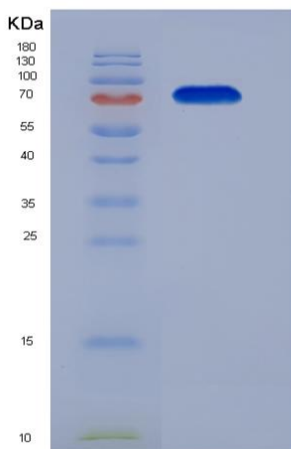


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.