

## Recombinant X-linked Inhibitor Of Apoptosis Protein (XIAP)

Catalog No.: TP07933

100µg

### Sequence Information

**Species:** Human

**Gene ID:**331

**Swiss Prot:**P98170

**Synonyms:**API3;BIRC4;hIAP-3;hIAP3;  
IAP-3;ILP1;MIHA;XLP2

**Residues:**Met1~Ser497

MTFNSFEGSKTCVPADINKEEEFVEEFNRLKTFANFSPGSPVSASTLARAGFLY  
TGEGDTVRCFSCHAAVDRWQYGD SAVGRHRKVSPNCRFINGFYLENSATQSTNS  
GIQNGQYKVENYLGSRDHFALDRPSETHADYLLRRTGQVVDISDTIYPRNPAMYS  
EEARLKS FQNWPDYAHLTPRELASAGLYYTGIGDQVQCFCCGGKLNWEPCDRA  
WSEHRRHFPCFFVLGRN LNIRSESDAVSSDRNFPNSTNLPRNPSMADYEARIF  
TFGTWIYSVNKEQLARAGFYALGEGDKVKCFHCGGGLTDWKPS EDPWEQHAKWY  
PGCKYLLEQKGQEYINNIHLTHSLEECLVRTTEKTPSLTRRIDDTIFQNP MVQE  
AIRMGFSFKDIKKIMEEKIQISGSNYKSLEVLVADLVNAQKDSMQDESSQTS LQ  
KEISTEEQLRRLQEELCKICMDRNIAIVFVPCGHLVTCKQCAEAVDKCPMCYT  
VITFKQKIFMS

### Product Information

**Source:** Recombinant expression.

**Host:** *E.coli*

**Tags:** N-terminal His-Tag

**Subcellular Location:** Nucleus, Cytoplasm.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.2

**Predicted Molecular Mass:** 60.4kDa

**Accurate Molecular Mass:** 60kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

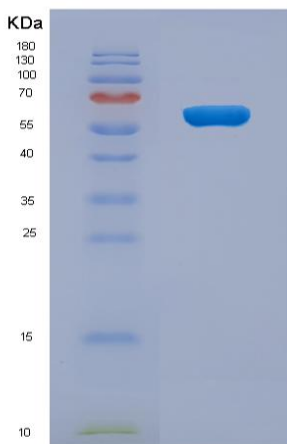


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.