



# Recombinant Peroxisome Proliferator Activated Receptor Gamma (PPAR $\gamma$ )

Catalog No.: TP09508

100 $\mu$ g

## Sequence Information

**Species:** Mouse

**Gene ID:**19016

**Swiss Prot:**P37238

**Synonyms:**PPAR-G; PPARG1; PPARG2; NR1C3

**Residues:**Met1-Tyr505

MGETLGDSPVDPEHGAFADALPMSTSQEITMVDTEMPFWPTNFGISSVDLSVMEDHSHSFDIKPFTTVDFSSISAPHYEDIPF  
TRADPMVADYKYDLKLQEYQSAIKVEPASPPYYSEKTQLYNRPHEEPSNSLMAIECRVCGDKASGFHYGVHACEGCKGFFRRT  
IRLKLIIYDRCDLNCRIHKKSRNKCQYCRFQKCLAVGMSHNAIRFGRMPQAEKEKLLAEISSDIDQLNPESADLRALAKHLYDS  
YIKSFPLTKAKARAILTGKTTDKSPFVIYDMNSLMMGEDKIKFKHITPLQEQSKEVAIRIFQGCQFRSVEAVQEITEYAKNIP  
GFINLDLNDQVTLTKYGVHEIIYTMLASLMNKDGVLISEGQGFMTREFLNLRKPFGDFMEPKFEFAVKFNALELDDSDLAIF  
IAVIIISGDRPGLLNKPIEDIQDNLLQALELQLKLNHPESSQLFAKVLQKMTDLRQIVTEHVQLLHVIKKTETDMSLHPLLQ  
EIIYKDLY

## Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** N-terminal His-Tag+Sumo

**Subcellular Location:** Secreted.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS (PH7.4) , containing 5% Trehalose.

**Original Concentration:** 200 $\mu$ g/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.9

**Predicted Molecular Mass:** 70.5kDa

**Accurate Molecular Mass:** 71kDa as determined by SDS-PAGE reducing conditions.

## [ USAGE ]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under



appropriate storage condition.

### [ IDENTIFICATION ]

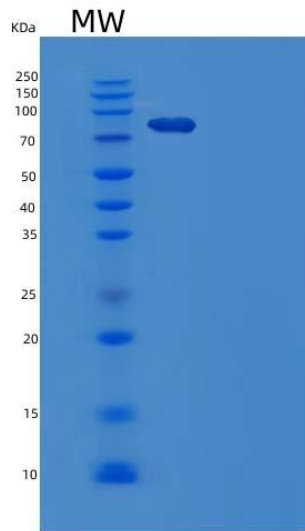


Figure . SDS-PAGE

### [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.