

Recombinant Mouse Interleukin 27A (IL27A)

Catalog No.: TP09623

50µg

Sequence Information

Species: Mouse

Gene ID:246779

Swiss Prot:Q8K3I6

Synonyms:p28; IL30; Interleukin 30;

Interleukin-27 subunit alpha; IL-27

subunit alpha; IL27-A

Residues:Phe29~Pro243

FPTDPLSLQELRREFTVSLYLARKLLSEVQGYVHSFAESRLPGVNLDLLPLGYH

LPNVSLTFQAWHHLSDSERLCFLATTLRPFPA MLGGLGTQGTWTS SEREQLWAM

RLDLRDLHRHLRFQVLAAGFKCSKEEDKEEEEEEEEEKKLPLGALGGPNQVS

SQVSWPQLLYTYQLLHSL ELVLSRAVRDLLLSLPRRPGSAWDS

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.6

Predicted Molecular Mass: 26.2Da

Accurate Molecular Mass: 27kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

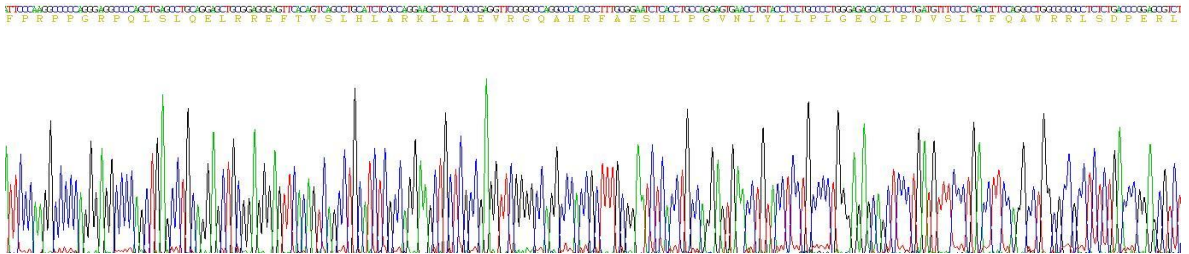


Figure 1. Gene Sequencing (Extract)

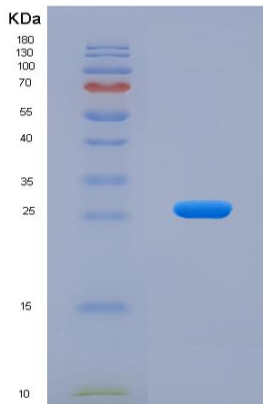


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.