

Recombinant Mouse Angiotensinogen (AGT)

Catalog No.: TP09561

50µg

Sequence Information

Species: Mouse

Gene ID:

Swiss Prot:P11859

Synonyms: ANHU; SERPINA8; Angiotensin;

Renin Substrate; Serpin Peptidase

Inhibitor, Clade A, Member 8 (Alpha-1

Antiproteinase, Antitrypsin)

Residues: Asp25~Val477

DRVYIHPFHLLYHNKSTCAQLENPSVETLPESTFEPVPIQAKTSPVNEKTLHDQ
LVLA AEKLEDEDRKRAAQVAMIANFVGFRMYKMLNEAGSGASGAILSPPALFGT
LVSFYLGSLDPTASQLQTLTDVPVKEGDCSTRLDGHKVLAAALRAVQGLLVLTQGG
SSSQTPLLQSIMVGLFTAPGFR LKHSFVQSLALFTPALFPRSLDLSTDPVLATE
KINRFIKAVTGWKMNLPLEGVSTDSTLLFNTRYVHFQGTMRGFSQLPGVHEFWVD
NSISVSVP MISGTGNFQHWSDAQNNFSVTCVPLGERATLLLIQPHCTSDLDRVE
ALIFRNDLLTWIENPPPRAIRLTLPQLEIRGSYNLQDLLAEDKLPTLLGAEANL
SNIGDTNPRVGEVLNSILLELKAGEEEQPTTSVQQPGSPEALDVTLSPPFLFAI
YEQDSGTLHFLGRVNNPQSVV

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.0

Predicted Molecular Mass: 51.2Da

Accurate Molecular Mass: 65kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

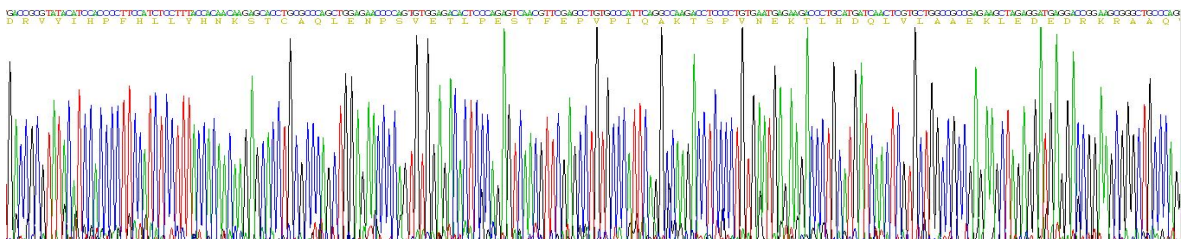


Figure 1. Gene Sequencing (Extract)

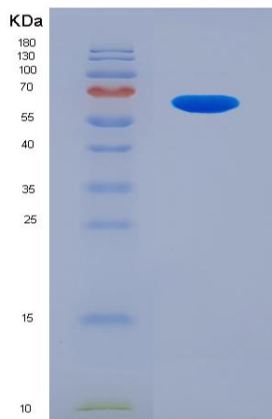


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

