

Recombinant Rat Interferon Alpha 11 (IFNa11)

Catalog No.: TP09565

50µg

Sequence Information

Species: Rat

Gene ID:

Swiss Prot:

Synonyms: Limitin

Residues: Thr20~Lys191

TCCLGCDLPHTHNLRNKRAFTLLAQMRRLSPVSLKDRQDFGFPLEKVVDG
QQIQKAQAIPVLHELTQQILSLFTSKESSTAWDATLLDSFCNDLQQQLSGLQA
CQMQQVGVQESHLAQEDSLLAVREYFHRITVYLREKKHSPCAWEVVRAEV
WRALSSSANLLARLREEKK

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: N-terminal His Tag

Subcellular Location: Secreted.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.2

Predicted Molecular Mass: 21.3kDa

Accurate Molecular Mass: 22kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was

30CTTCGTCAGATGGCCTCCACACTCAATCGAAGAGGCATTGCCCTCCAAATGGGCTCCCCTGTCTTGGGAGAGGAGTTTGATATCTTGGAAGCATGAGCATGAGAGCTCAATCCGCTCTGGCTGGAGAGAATCAATCCCTCTCATAGGGCACTACTCTGGA
 TCCLGCDLPHTNLRNKRAEFTLLAQMRRILSPVSCCLKDRDEFGFPLEKVDGQOTOKAQAITPVLTHELTOOITLSLETSSKESSSTAMU

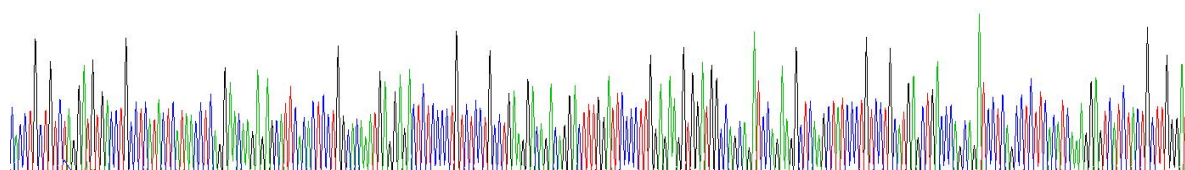


Figure 1. Gene Sequencing (Extract)

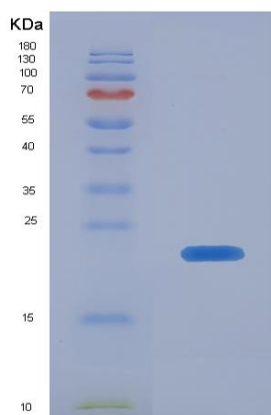


Figure 2. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.