

Recombinant Human A Disintegrin And Metalloprotease 8 (ADAM8)

Catalog No.: TP09646

50µg

Sequence Information

Species: Human

Gene ID:101

Swiss Prot:P78325

Synonyms:CD156-A; MS2; CD156a;

Disintegrin and metalloproteinase
domain-containing protein 8

Residues:Met1-Pro497

MRGLGLWLLGAMMLPAIAPSRPWALMEQYEVVLPWRLPGPRVRRALPSHLGLHP
ERVSIVLGATGHNFTLHLRKNRDLLGSGYTETYTAANGSEVTEQPRGQDHCIFYQ
GHVEGYPDSAASLSTCAGLRGFFQVGSDDLHLEPLDEGGEGGRHAVYQAEHLLQ
TAGTCGVSDDSLGLGPRTAAVFRPRPGDSLPSRETRYVELYVVVDNAEFQML
GSEAAVRHRVLEVVNHVDKLYQKLNFRVVLVGLFIWNSQDRFHVSPDPSVTLEN
LLTWQARQRTRRHLHDNVQLITGVDFGTGTVGFARVSAMCSHSSGAVNQDHSKN
PVGVACTMAHEMHGNLGMHDENVQGCRQCERFEAGRCIMAGSIGSSFPRMFSD
CSQAYLESFLERPQSVCLANAPDLSHLVGGPVCGNLFVERGEQCDGPPEDCRN
RCCNSTTCQLAEGAQCAHGTCCECKVKPAGELCRPKKDMCDLEEFCDGRHPEC
PEDAFQENGTP

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.1

Predicted Molecular Mass: 60kDa

Accurate Molecular Mass: 60kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

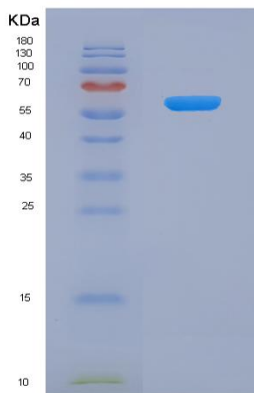


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.