

Recombinant Human Glial Cell Line Derived Neurotrophic Factor Receptor Alpha 2 (GFRA2)

Catalog No.: TP04611

50µg

Sequence Information

Species: Human

Gene ID:2675

Swiss Prot:O00451

Synonyms: GDNF Family Receptor Alpha-2;

GDNF Receptor Alpha-2;

GDNFR-Alpha-2; GFR-Alpha-2; GDNF

Receptor Beta; GDNFR-Beta; Neurturin

Receptor Alpha; NRTNR-Alpha;

NTNR-Alpha; RET Ligand 2;

TGF-Beta-Related Neurotrophic Factor

Receptor 2; GFRA2; GDNFRB; RETL2;

TRNR2

Residues:Ser22-Ser441

SPSSLQGPPELHGWRPPVDCVRANELCAAESNCSSRYRTL RQCLAGRDRNTMLAN
KECQAALVLEVLQESPLYDCRCKRGMKKELQCLQIYWSIHLGLTEGEEFYEASPYE
PVTSRLSDIFRLASIFSGTGADPVVSAKSNHCLDAAKACNLNDNCKKLRSSYIS
ICNREISPTERCNRRKCHKALRQFFDRVPSEYTYRMLFCSCQDQACAERRRQTI
LPSCSYEDKEKPNCLDLRGVCRDHLCSRLADFHANCASYQTVTSCPADNYQ
ACLSYAGMIGFDMTPNYVDSSPTGIVSPWCSCRGSNMEECEKFLRDFTEN
PCLRNAIQAFNGTGVNVPKGPSFQATQAPRVEKTPSLPDDLSDSTSLGTSVI
TTCTSVQEQGLKANNSKELSMCFTELTTNIIPGSNKVIKPN

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.5

Predicted Molecular Mass: 70kDa

Accurate Molecular Mass: 70kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

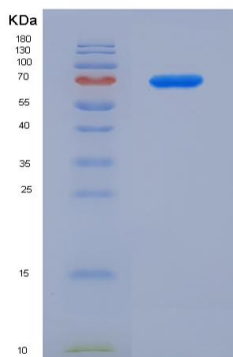


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.