

Recombinant Human Interleukin 11 Receptor Alpha (IL11Ra)

Catalog No.: TP08244

50µg

Sequence Information

Species: Human

Gene ID:3590

Swiss Prot:Q14626

Synonyms:CRSDA;hCG_2011440;IL11

RA;MGC2146

Residues:Ser24-Val363

SPCPQAWGPPGVQYGQPGRSVKLCCPGVTAGDPVSWFRDGEPKLLQGPDSSGLGH

ELVLAQADSTDEGTICQTLDGALGGTVTLQLGYPPARPVVSCQAADYENFSCT

WSPSQISGLPTRYLTYSRKKTVLGADSQRRSPSTGPWPCPDPLGAARCVVHGA

EFWSQYRINVTEVNPLGASTRLLDVSLQSILRPDPPQGLRVESVPGYPRRLRAS

WTYPASWPCQPHFLKFRLLQYRPAQHPAWSTVEPAGLEEVITDAVAGLPHAVRV

SARDFLDAGTWSTWSPEAWGTPSTGTIPKEIPAWGQLHTQPEVEPQVDSPPPR

PSLQPHPRLLDHRDSV

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Cell membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.5

Predicted Molecular Mass: 40kDa

Accurate Molecular Mass: 40kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

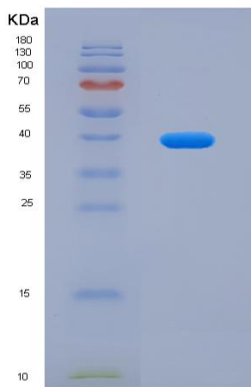


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.