

Recombinant Human Insulin Like Growth Factor 2 Receptor (IGF2R)

Catalog No.: TP09597

50µg

Sequence Information

Species: Human

Gene ID:3482

Swiss Prot:P11717

Synonyms:CD222; IGF2-R; CIMPR; M6P-R;

MPRI; Cation-Independent Mannose-6

Phosphate Receptor; M6P/IGF2 receptor;

300 kDa mannose 6-phosphate receptor

Residues:Ala43-Gly1365

APFPELCSYTWEAVDTKNNVLYKINICGSVDIVQCGPSSAVCMHDLKTRTYHSV
GDSVLRSATRSLLEFNTTVSCDQQGTNHRVQSSIAFLCGKTLGTPEFVTATECV
HYFEWRTTAACKKDIFKANKEVPCYVFDEELRKHDNLPLIKLSGAYLVDDSDPD
TSLFINVCRDIDTLRDPGSQLRACPPGTAACLVRGHQAFDVGQPRDGLKLVRKD
RLVLSYVREEAGKLDFCDGHSPAVTITFVCPSEERREGTIPKLTAKSNCRYEIEW
ITEYACHRDYLESKTCSLSGEQQDVSIDLTPLAQSGGSSYISDGKEYLFYLNVC
GETEIQFCNKKQAAVCQVKKSDTSQVKAAGRYHNQTLRYSDGDLTLIYFGGDEC
SSGFQRMSSVINFECKNTAGNDGKGTVPFTGEVDCTYFFTWDEYACVKEKEDLL
CGATDGKKRYDLSALVRHAEPEQNWEAVDGSQTETEEKHFFINICHRVLQEGKA
RGCPEDAAVCAVDKNGSKNLGKFISSPMKEKGNIQLSYSDGDDCGHGKKIKTNI
TLVCKPGDLESAPVLRSTGEGGCFYEFEWHTAAACVLSKTEGENCTVFDSQAGF
SFDLSPLTKKNGAYKVETKKYDFYINVCGPVSVSPCQPDGACQVAKSDEKTWN
LGLSNAKLSYDGMIIQLNRYGGTPYNNERHTPRATLITFLCDRDAGVGFPEYQE
EDNSTYNFRWYTSYACPEEPLCCVVDTPSTLEQYDLSSLAKSEGGGLGWNWYAMD
NSGEHVTWRKYVINVCRLNPVPGCNRYASACQMKYEKDGQSFTVVVISNLGM
AKTGPVVEDSGSLLLEYVNGSACTTSDGRQTTYTTRIHLVCSRGRNLNSHPIFSL
NWEVVSFLWNTAEACPIQTTTDDQACSIRDPNSGFVFNLNPLNSSQGYNVSG
IGKIFMFNVCGTMPVCGTILGKPASGCEAETQTEELKNWKPAPVVGIEKSLQLS
TEGFITLTYKGPLSAKGTADAFIVRFVCNDDVYSGPLKFLHQDIDSGQGIRNTY
FEFETALACVPSPVDCQVTDLAGNEYDLTGLSTVRKPWTAVDTSVDGRKRTFYL
SVCNPLPYIPGCQGSAVGSCLVSEGNLWNLGVVQMSPQAAANG

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Secreted.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 5.1

Predicted Molecular Mass: 160kDa

Accurate Molecular Mass: 160kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

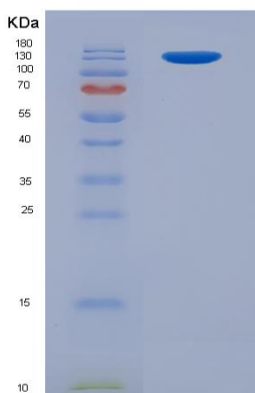




Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.