



Recombinant Human Tumor Necrosis Factor Receptor Superfamily, Member 21 (TNFRSF21)

Catalog No.: **TP04612** **50µg**

Sequence Information

Species: Human

Gene ID:27242

Swiss Prot:O75509

Synonyms:Tumor Necrosis Factor Receptor

Superfamily Member 21; Death Receptor

6; CD358; TNFRSF21; DR6

Residues:Met1-Leu350

MGTSPSSSTALASCSRIARRATATMIAGSLLLLGFLSTTTAQPEQKASNLIGTY
RHVDRATGQVLTCDKCPAGTYVSEHCTNTSLRVCSSCPVGTFTRHENGIEKCHD
CSQPCPWPMIEKLPAAALTDRECTCPGMEFQSNATCAPHTVCPVGWGVRRKKGTE
TEDVRCKQCARGTFSDVPSSVMKCKAYTDCLSQLLVVIKPGTKETDNVCGTLPS
FSSSTSPSPGTAIFRPEHMETHEVPSSTYVPKGMNSTESNSSASVRPKVLSSI
QEGTVPDNTSSARGKEDVNKTLPNLQVVNHQQGPHHRHILKLLPSMEATGGEKS
STPIKGPGRGHPRQNLHKHFDINEHL

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag and Fc

Subcellular Location: Membrane.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.3

Predicted Molecular Mass: 95kDa

Accurate Molecular Mass: 95kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

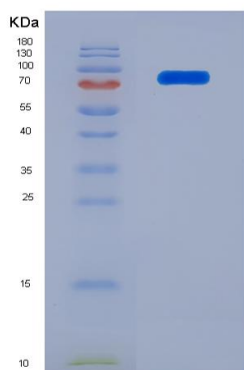


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.