

Recombinant Human Galactosidase Beta (GLb)

Catalog No.: TP04609

50µg

Sequence Information

Species: Human

Gene ID:2720

Swiss Prot:P16278

Synonyms:Beta-Galactosidase; Acid

Beta-Galactosidase; Lactase; Elastin

Receptor 1; GLB1; ELNR1

Residues:Met1-Val677

MPGFLVRILPLLLVLLLLGPTRGLRNATQRMFEIDYSRDSFLKDGQPFRYISGS
IHYSRVPRFYWKDRLLKMKMAGLNAIQTYVPWNFHEPWPGQYQFSEDHDEVEYFL
RLAHELGLLVILRPGPYICAEWEMGGLPAWLLEKESILLRSSDPDYLAADV KWL
GVLLPKMKPLLYQNGGPVITVQVENEYGSYFACDFDYLRFLQKRFRHHLGDDVV
LFTTDGAHKTF LKCGALQGLYTTVDFTGSGNITDAFLSQRKCEPKGPLINSEFY
TGWLDHWGQPHSTIKTEAVASSLYDILARGASVNLYMFIGGTNFAYWNGANSPY
AAQPTS YDYDAPLSEAGDLTEKYFALRNIIQKF EKVPEGPIPPSTPKFAYGKVT
LEKLKTVGAALDILCPSPGIKSLYPLTFIQVKQHYGFVLYRTTLPQDCSNPAPL
SSPLNGVHDRAYVAVDGIPQGVLERNNVITL NITGKAGATLDLLVENMGRVNYG
AYINDFKGLVSNLT LSSNILT DWTIFPLDTEDAVRSHLGGWGHRS GHHDEAWA
HNSSNYTLPAFYMGNFSIPSGIPDL PQDTFIQFPGWTKGQVWINGFNLGRYWPA
RGPQLTLFVPQHILMTSAPNTITVLELEWAPCSSDDPELCAVTFVDRPVIGSSV
TYDHPSKPVEKRLMPPPPQKNKDSWLDHV

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Lysosome .

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.6

Predicted Molecular Mass: 74.6kDa

Accurate Molecular Mass: 75kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

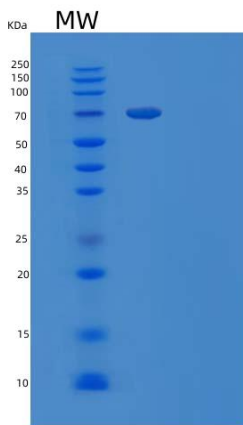


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.