

Recombinant Human Hepatocyte Growth Factor (HGF)

Catalog No.: TP09698

50µg

Sequence Information

Species: Human

Gene ID:15234

Swiss Prot:Q08048

Synonyms:HGF; F-TCF; HGFB; HPTA; SF;

Scatter Factor; Hepapoinetin A;

Fibroblast-Derived Tumor Cytotoxic

Factor; Lung Fibroblast-Derived

Mitogen

Residues:Gln33-Lys728

QKKRRNTLHEFKKSAKTTLTKE DPLLKIKTKKVNSADECANRCIRNRGFTFTCK
AFVFDKSRKRCYWYPFNSMSSGVKKGFGHEFDLYENKDYIRNCIIIGKGGSYKGT
VSITKSGIKCQPWNSMIPHEHSFLPSSYRGKDLQENYCRNPRGEEGPGWCFTSN
PEVRYEVC DIPCSEVEECMTNGESYRGPM DHTESGKTCQRWDQQT PHRHKFLP
ERYPDKGFD DNYCRNPDGKPRPWCYTLD PDPWEYCAIKTCAHS AVNETDVPME
TTECIQGQGEGRGTSNTI WNGIPCQRWDSQYPHKHDITPENFKCKDLRENYCR
NPDGAESPWCFTTDPNIRVGYCSQIPKCDVSSGQDCYRGN GKNYMGNLSKTRSG
LTCSMWDKNMEDLHRHIFWEPDASKLNKNYCRNPDDDAHGPWCYTGNPLIPWDY
CPISRCEGDTTPTIVNL DHPVISCAKTKQLRVVNGIPTQT TVGWMVSLKYRNKH
ICGGS LIKESWLTARQCFPARNKDLKDYEAWLGIHDVHERGEEKRKQILNISQ
LVYGPESDLVLLKLARPAILDNFVSTIDLPSYGCTIPEKTTCSIYGWGYTGLI
NADGLLRVAHLYIMGNEKCSQHHQGVTLNESEL CAGAEKIGSGPCEGDYGGPL
ICEQH KM RMVLGVIVPGRGCAIPNRPGIFVRVAYYAKW IHKVILTYKL

Product Information

Source: Eukaryotic expression.

Host: 293F cell

Tags: C-terminal His-Tag

Subcellular Location: Extracellular matrix.

Purity: >95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

Original Concentration: 1000µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.0

Predicted Molecular Mass: 90kDa

Accurate Molecular Mass: 90kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

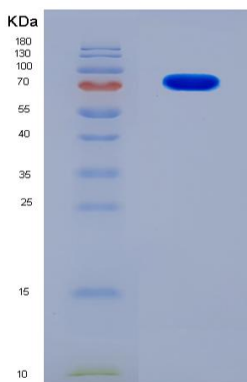


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.

