



Recombinant Plasminogen Activator Inhibitor 1 (PAI1)

Catalog No.: TP09740 200µg

Sequence Information

Species: Rat

Gene ID:24617

Swiss Prot:P20961

Synonyms:SERPINE1; PLANH1

Residues:Ser24-Pro402

SPLPESHTAQQATNFGVKVFQHVQASKDRNVVFSYGVSSVLAMLQLTTAGKTRQQIQDAMGFNISERGTAPALRKLSKELMGSWNKNEISTA
DAIFVQRDLELVQGFMPHFVKLFRTTVKQVDFSEVERARFIINDWVERHTKGMISDLLAKGAVNELTRLVLVNALYFNGQWKTPFLEASTHQRL
FHKSDGSTISVPMMAQNNKFNYTEFTTPDGHEYDILELPYHGETLSMFIAAPFEKDVPLSAITNILDALIRQWKSNMTRLPRLLILPKFSLET
EVDLRGPLEKLGMTDIFSSTQADFTSLSDQEQLSVAQALQVKIEVNESGTVASSSTAILVSARMAPTEMVLDRSFLFVVRHNPTETILFMGQL
MEP

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His

Subcellular Location: Membrane.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.2, containing 0.01% skl,5%Trehalose.

Original Concentration: 100µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.5

Predicted Molecular Mass: 46kDa

Accurate Molecular Mass: 46kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

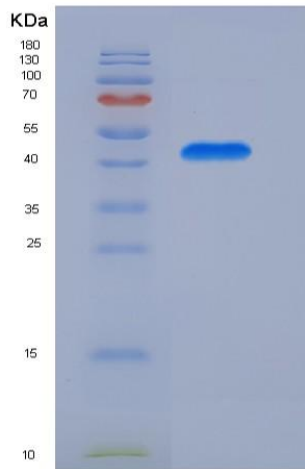


Figure 1. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.