

## Recombinant Mouse Interleukin-12 Subunit $\beta$ (IL-12B)

Catalog No.: TP05328

50 $\mu$ g

### Sequence Information

**Species:** Mouse

**Gene ID:**16160

**Swiss Prot:**P43432

**Synonyms:**Interleukin-12 subunit beta; IL-12B;  
Cytotoxic lymphocyte maturation factor  
40 kDa subunit; CLMF p40; IL-12  
subunit p40; Il12b

**Residues:**Met23-Ser335

MWELEKDVYVVEVDWTPDAPGETVNLTCDTPEEDDITWTSQQRHGVIGSGKTLT  
ITVKEFLDAGQYTCHKGGETLSHSHLLLHKKENGIWSTEILKNFKNKTLKCEA  
PNYSGRFTCSWLVRNMDLKFNIKSSSSSPDSRAVTCGMASLSAEKVTLQDRDY  
EKYSVSCQEDVTCPTAEETLPIELALEARQQNKYENYSTSFFIRDIIKPDPPKN  
LQMKPLKNSQVEVSWEYPSWSTPHSYFSLKFFVRIQRKKEKMKETEEGCNQKG  
AFLVEKTSTEVQCKGGNVCVQAQDRYNSSCSKWACVPCRVR

### Product Information

**Source:** Eukaryotic expression.

**Host:** 293F cell

**Tags:** C-terminal His Tag

**Subcellular Location:** Secreted.

**Purity:** >95%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and Proclin300.

**Original Concentration:** 1000 $\mu$ g/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.3

**Predicted Molecular Mass:** 47.1kDa

**Accurate Molecular Mass:** 47kDa as determined by SDS-PAGE reducing conditions.

### [ USAGE ]

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-0.5 mg/mL. Do not vortex.

## [ STORAGE AND STABILITY ]

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

## [ IDENTIFICATION ]

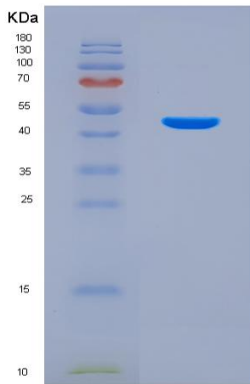


Figure 1. SDS-PAGE

## [ IMPORTANT NOTE ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.