

Recombinant Neurofilament, Light Polypeptide (NEFL)

Catalog No.: TP09905

50µg

Sequence Information

Species: Human

Gene ID:4747

Swiss Prot:P07196

Synonyms: CMT1F; CMT2E; NF-L; NF68; NFL; 68 kDa neurofilament protein; Neurofilament triplet L protein

Residues:Ser2-Asp543

SSFSYEPYYSTSYKRRYVETPRVHISSVRSGYSTARSAYSSYSAPVSSSLSVRRSYSSSSGSLMPSELENLDLSQVAAISND
LKSIRTQEKAQLQDLNDRFASFIERVHELEQQNKVLEAELLVLRQKHSEPSRFRALYEQEIRDLRLAAEDATNEKQALQGE
REGLEETLRNLQARYEEEVLSREDAEGRLEARKGADEAALARAELEKRIDSLMDEISFLKKVHEEEIAELQAQIQYAQIS
VEMDVTKPDLSAALKDIRAQYEKLAAKNMQNAEEWFKSRFTVLTESAAKNTDAVRAAKDEVSESRLLKAKTLEIEACRGM
NEALEKQLQELEDKQNADISAMQDTINKLENELRTTKSEMARYLKEYQDLLNVKMALDIEIAAYRKLLGEETRLSFTSVG
SITSGYSQSSQVFGRSAYGGLQTSSYLMSTRSFPSSYYTSHVQEEQIEVEETIEAAKAEAAKDEPPSEGEAEKKDEKKEAE
EEEEAEFFEEAAKEESEEAKEEFEGGEGEGEGEETKEAEKKKVEGAGEEOAAKKKD

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-terminal His-Tag

Subcellular Location: Secreted.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS (PH7.4) , containing 5% Trehalose.

Original Concentration: 400µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.3

Predicted Molecular Mass: 65kDa

Accurate Molecular Mass: 70kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.



[IDENTIFICATION]

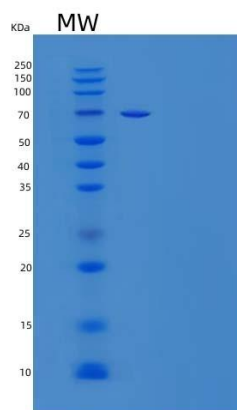


Figure . SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.