



## Recombinant Human CCNB1 Protein

Catalog No.: TP01472

1mg

### Sequence Information

**Species:** Human

**Gene ID:**891

**Swiss Prot:**P14635

**Synonyms:**Cyclin B1, CCNB

**Residues:**Met1-Val433

MALRVTRNSKINAENKAKINMAGAKRVPTAPAATSKPGLRPRALTALGDIGNKVSEQLQAKMPMKKEAKPSATGKVIDKKLPKPLEKVPMLVPVPVSEP  
VPEPEPEPEPEPVKEEKLSPPEILVDTASPSMETSGCAPAEEDLCQAFSDVILAVNDVDAEDGADPNLCSEYVKDIYAYLRQLEEEQAVRPKYLLG  
REVTGNMRAILIDWLQVQMKFRLLQETMYMTVSIIDRFMQNNCPVKMMLQLVGVTAMFIASKYEEMYPPEIGDFAFVTDNTYTKHQIRQMEMKILR  
ALNFGGLGRPLPLHFLRRASKIGEVDVEQHTLAKYLMELTMDYDMVHFPPSQIAAGAFCLALKILDNGEWPTLQHYLSYTEESLLPVMQHLAKNVV  
MVNQGLTKHMTVKNKYATSKHAKISTLPQLNSALVQDLAKAVAKV

### Product Information

**Source:** Prokaryotic expression.

**Host:** *E. coli*

**Tags:** N-6His

**Subcellular Location:** Nucleus, Cytoplasm, Chromosome.

**Purity:** >90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, 5%Trehalose .

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 8.02

**Predicted Molecular Mass:** 51.1kDa

**Accurate Molecular Mass:** 51kDa as determined by SDS-PAGE reducing conditions.

### **[ USAGE ]**

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

### **[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

### **[ IDENTIFICATION ]**

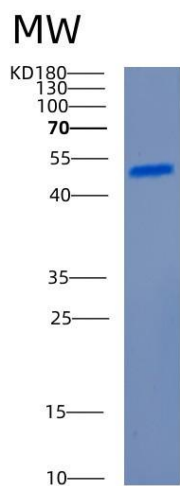


Figure . SDS-PAGE

### [ **IMPORTANT NOTE** ]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.