



Recombinant Human CCNB1 Protein

Catalog No.: TP09892

1mg

Sequence Information

Species: Human

Gene ID:891

Swiss Prot:P14635

Synonyms:Cyclin B1, CCNB

Residues:Met1-Val433

MALRVTRNSKINAENKAKINMAGAKRVPTAPAATSKPGLRPR TALGDIGNKVSEQLQAKMPMKKEAKPSATGKVIDKKLPKPLEKVPMLVPVPVSEP
VPEPEPEPEPEPVKEEKLSP EPI LVDTASPSMETSGCAPAEEDLCQAFSDVILAVNDVDAEDGADPNLCSEYVKDIYAYLRQLEEEQAVRPKYLLG
REVTGNMRAILIDWL VQVQMKFRLLQETMYMTVSIIDRFMQNNCVPKKMLQLVGVTAMFIASKYEEMYPPEIGDFAFVTDNTYTKHQIRQMEMKILR
ALNFG LGRPLPLHFLRRASKIGEVDVEQHTLAKYLMELTMDYDMVHFPPSQIAAGAFCLALKILDNGEWTP TLQHLYSYTEESLLPVMQHLAKNVV
MVNQGLTKHMTVKNKYATSKHAKISTLPQLNSALVQDLAKAVAKV

Product Information

Source: Prokaryotic expression.

Host: *E. coli*

Tags: N-6His+GST

Subcellular Location: Nucleus, Cytoplasm, Chromosome.

Purity: >90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, 5%Trehalose .

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.02

Predicted Molecular Mass: 77.1kDa

Accurate Molecular Mass: 77kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in PBS (pH7.4) to a concentration of 0.1-0.5 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]

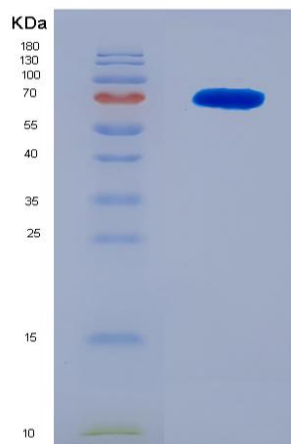


Figure . SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for in vitro and research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.